




Gorilla® Evans Osteotomy

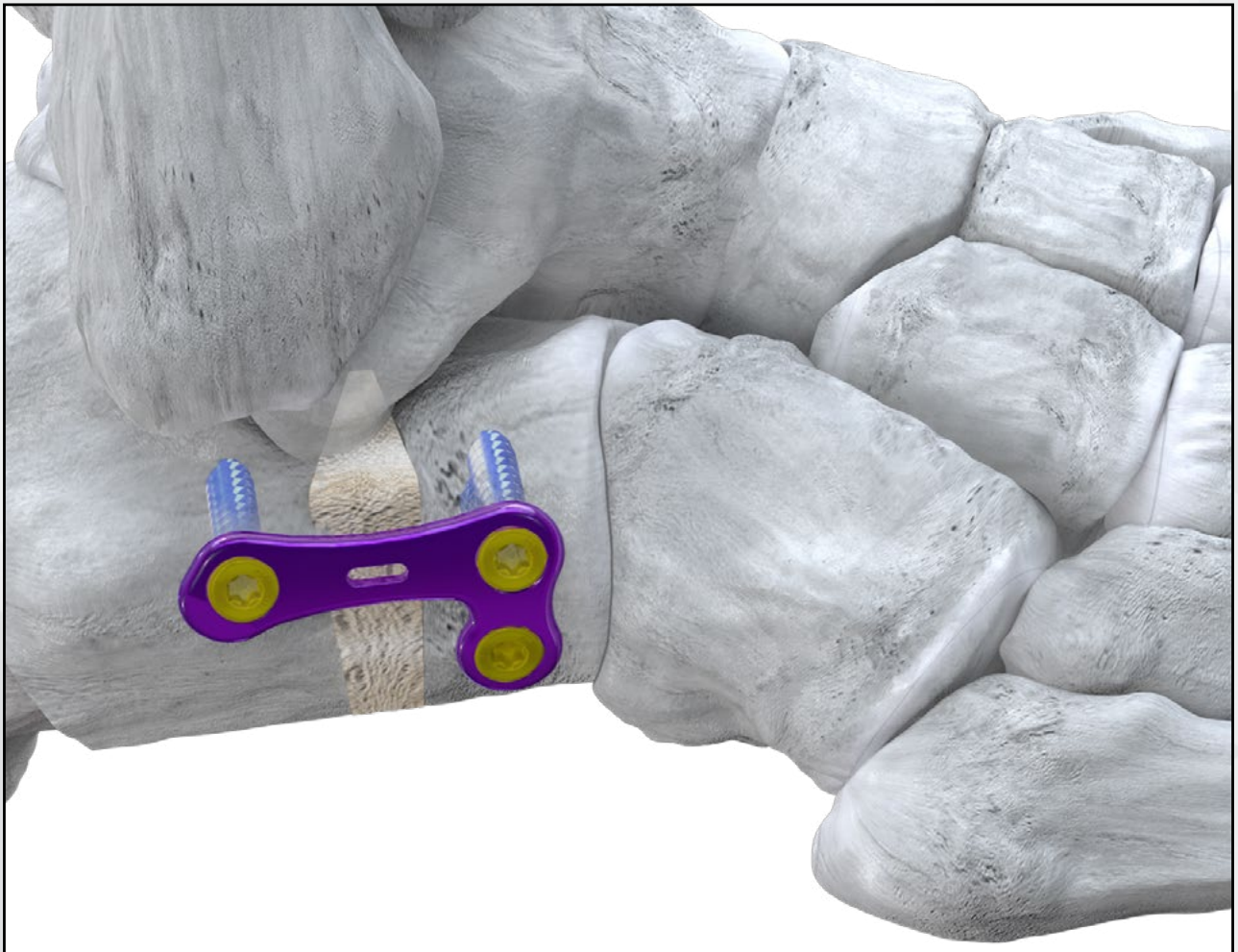

PRESERVE™ EVANS LATERAL COLUMN LENGTHENING GRAFT

Patented Wedge Geometry

- ▶ Dorsal to plantar taper helps reduce risk of attenuation of the lateral aspect of the long plantar ligament
- ▶ Lateral to medial taper helps fit the shape of the Evans osteotomy and allows for easier insertion

Aseptically Processed

- ▶ Processed without gamma irradiation or hydrogen peroxide to help preserve the native mechanical advantages of human bone and the osteoinductivity of the environment in which the graft is being implanted

Density Matched

- ▶ Donor harvest site matched specific to Evans procedure for optimal strength and blood flow requirements
 - Donor sites include: Patella, talus and femoral calcar



PIN WIRE DISTRACTOR

- ▶ Sized appropriately for foot and ankle applications
- ▶ Can create up to a 28 mm opening
- ▶ Ratcheting handle locks to maintain rigid distraction across surgical site
- ▶ Included in every Gorilla® Plating System



HEVANS® PLATE

Low Profile


- 1.1 mm plate design

Peroneal Tendon Ramp

- Ramped proximal arm designed specifically to help reduce irritation of the peroneus brevis and surrounding soft tissue

Stable Fixation

- Design and orientation of the plate intended to stabilize anterior calcaneal fragment and maintain correction during graft incorporation

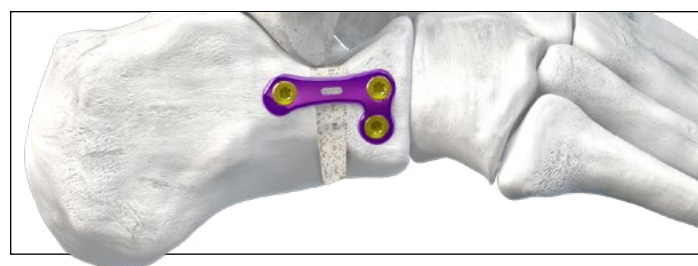

HEvans® plate used in conjunction with PRESERVE™ Evans Graft

BOW & ARROW® EVANS WEDGE PLATE

4 symmetrical size options for use in left or right foot




Plantar soft tissue ramps help to reduce irritation to peroneal tendons

Arrow on "bow" hooks around proximal cortex to help prevent lateral dislocation of plate



Tapered "bow" helps ease insertion and offload surrounding soft tissue structures

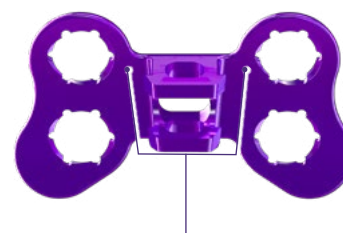

Built-in osteotomy spacer matches the geometry of PRESERVE™ Evans Graft



	HEVANS® PLATE	EVANS BOW & ARROW® PLATE
		
Total Plates	6	4
Plate Thickness	1.1 mm	1.1 mm
Size	Small, Medium, Large	6 mm, 8 mm, 10 mm, 12 mm (Spacer Size)
Left and Right Offerings	Yes	Symmetrical



Gorilla® Evans Osteotomy

EVAN-01 RevD
2023-10-05

Paragon 28®, Inc. 
14445 Grasslands Dr.
Englewood, CO 80112 USA
(855) 786-2828

Paragon 28® Medical Devices Trading Limited
First Floor Block 7 Beckett Way
Park West Business Park
Dublin 12
D12 X884
Ireland
+353 (0) 1588 0350

™Trademarks and ®Registered Marks of Paragon 28®, Inc.
© Copyright 2023 Paragon 28®, Inc. All rights reserved.
Patents: www.paragon28.com/patents

Exclusively foot & ankle ²⁸
Paragon®

www.Paragon28.com

For the contraindications, potential complications and adverse reactions, warnings and precautions associated with this device, please refer to the device specific instructions for use at <http://www.paragon28.com/ifus>