




**Gorilla<sup>®</sup> R3LEASE<sup>™</sup> Stabilization System**

**PRODUCT INFORMATION**

**Designed to Provide a Simple Solution for Syndesmosis Repair**

- Once a patient begins weightbearing, traditional syndesmotic screws may fracture.
- The location of screw fracture of a traditional syndesmotic screw can vary, leading to difficult removal situations and/or potential pain at the screw fracture site.
- The Gorilla<sup>®</sup> R3LEASE<sup>™</sup> Stabilization System was designed so that if screw fracture occurs, the screw breaks cleanly at the notch point in the clear space.

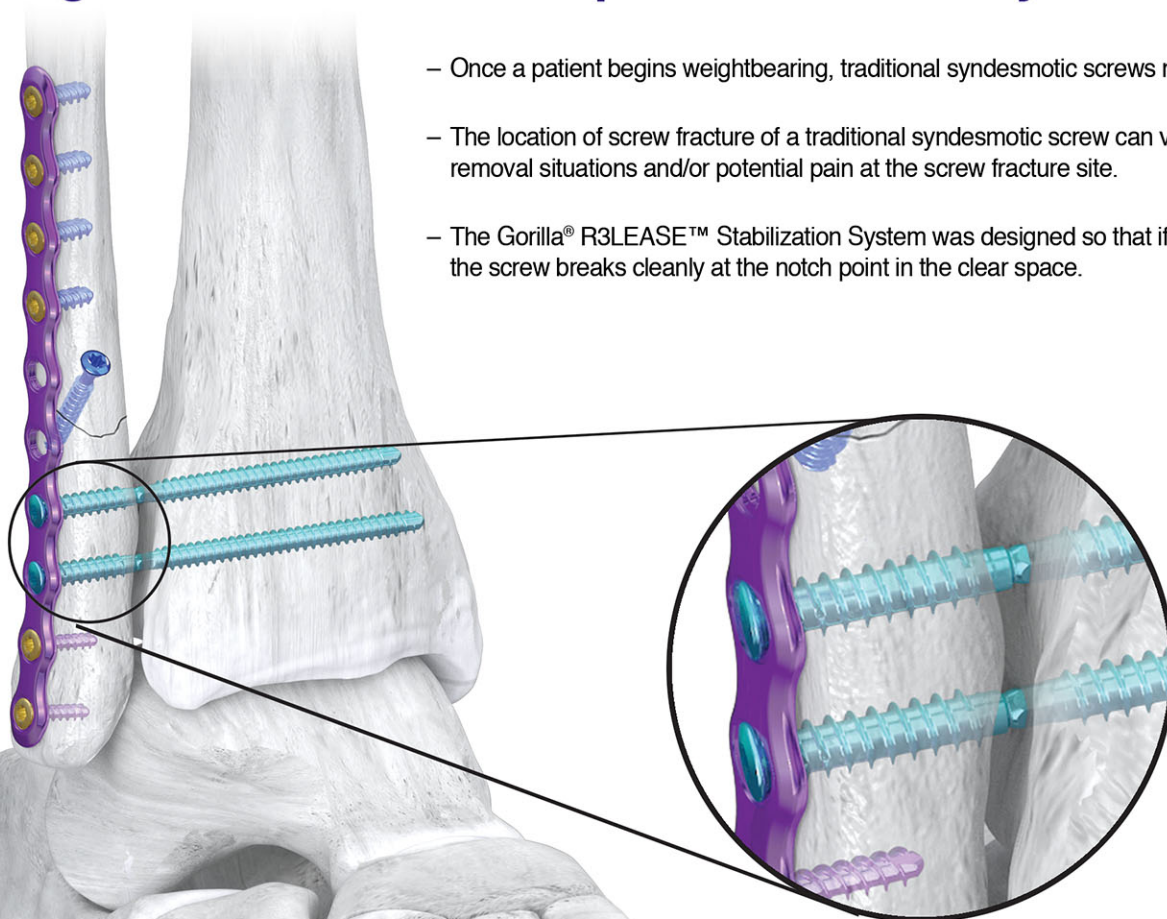

R3LEASE<sup>™</sup> notch resides in the clear space between fibula and tibia

# Gorilla® R3LEASE™ Stabilization System

## IMPLANT OFFERING




### R3LEASE™ SCREW SIZES

3.9 mm R3LEASE™ Screw - 14 mm Notch Length



3.9 mm R3LEASE™ Screw - 17 mm Notch Length

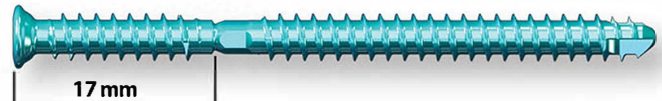


**Screw Lengths:**  
 40 mm - 45 mm (in 5 mm increments)  
 48 mm - 64 mm (in 2 mm increments)

**Screw Lengths:**  
 40 mm - 45 mm (in 5 mm increments)  
 48 mm - 64 mm (in 2 mm increments)

### SCREW FEATURES

- **Multiple Notch Lengths: 14 mm and 17 mm**
  - Addresses differences in patient anatomy and accommodates use with or without a plate
  - 2 mm increments from 48 - 64 mm intended for quadricortical fixation
- **Strength in Design**
  - Static bending strength and stiffness of the Gorilla® R3LEASE™ Screw exceeded that of a comparable 3.5 mm solid screw as well as a 4.0 mm cannulated screw made of the same material.<sup>1</sup>
- **If removal is desired, multiple removal features have been incorporated:**

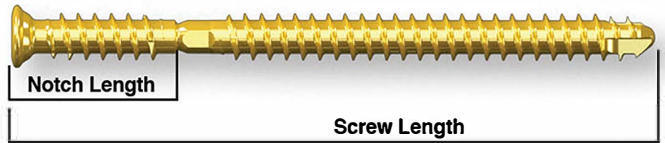

- **Lateral removal feature**
  - Hex feature allows for lateral removal of tibial fragment through the fibula

- **Medial removal feature**
  - Medial removal feature on screw tip allows for medial removal using a specialized driver

1. Internal data on file, TR-17102601

GBSS-01 RevD  
 ™ Trademarks and ® Registered Trademarks of Paragon 28, Inc.  
 © Copyright 2020 Paragon 28, Inc. All rights reserved.  
 Patents: [www.paragon28.com/index.php/patents](http://www.paragon28.com/index.php/patents)



- **Washer**
  - Available when the R3LEASE™ screw is used without or outside of a plate

Paragon 28, Inc.   
 14445 Grasslands Dr.  
 Englewood, CO 80112 USA  
 (855) 786-2828

Paragon 28 Medical Devices Trading Limited  
 First Floor Block 7 Beckett Way  
 Park West Business Park  
 Dublin 12  
 D12 X884  
 Ireland  
 +353 (0) 1588 0350

For additional information on Paragon 28® and its products please visit: <http://www.paragon28.com>

