



# Monster® BITE Screw System

## PRODUCT INFORMATION

#### **DESCRIPTION**

To address varying patient anatomies and address indications throughout the foot and ankle, the Monster® BITE Screw system includes 22 total implants in 2.0 mm and 2.7 mm diameters constructed from Titanium Alloy Ti 6AI-4V.

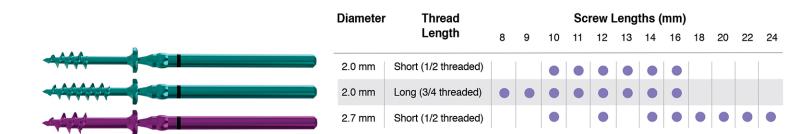

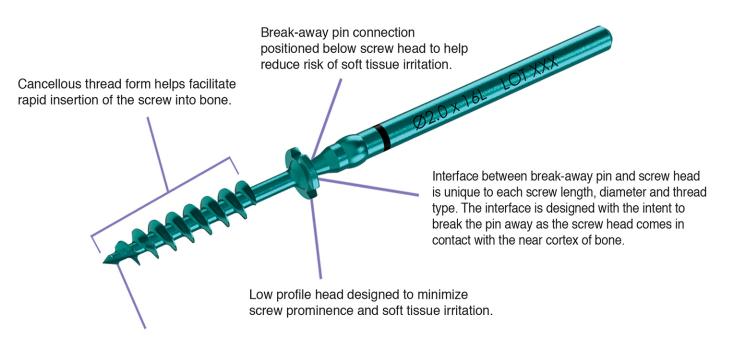

Self-drilling and self-tapping screw tip allows for screw insertion without a pilot hole.





# Monster® BITE Screw System

### SYSTEM SPECIFIC INFORMATION



### Monster® BITE 3 Prong Driver Shaft, A/O Connection

Aids in fully seating the screw in cases where the pin breaks away from the head of the screw prior to the screw head becoming flush to the bone.

Note: Each implant has a unique torque break away value calibrated for the length, diameter and thread length of the screw to help provide optimal break-away timing of the pin.



#### Weil Clamp

The Weil Clamp is designed to be placed around the osteotomy to temporarily secure the position of the distal fragment. The Weil Clamp can determine the approximate height of the metatarsal head by reading the distance off of the ratchet on the clamp.

#### MBSS-01 RevD

© Copyright 2020 Paragon 28, Inc. All rights reserved. Patents: www.paragon28.com/index.php/patents

Trademarks and Registered Trademarks of Paragon 28, Inc.

Paragon 28 Medical Devices Trading Limited First Floor Block 7 Beckett Way Park West Business Park Dublin 12 D12 X884 Ireland +353 (0) 1588 0350



For the contraindications, potential complications and adverse reactions, warnings and precautions associated with this device, please refer to the device specific instructions for use at http://www.paragon28.com/ifus

Paragon 28, Inc. 14445 Grasslands Dr. Englewood, CO 80112 USA (855) 786-2828



