

**Gorilla<sup>®</sup> MTP Plating System**

**PRODUCT INFORMATION**

**PLATING SYSTEM FEATURES**

- 32 plating options available in Short, Primary, Revision and Graft Spanning options
- Available in 0, 5, and 10 degrees of dorsiflexion (see right)
- Plates are 1.3 – 1.6 mm thick and machine contoured Ti-6Al-4V ELI
- Tightened distal cluster of screws to better match anatomy of the proximal phalanx
- Added two-hole distal plate for particularly small proximal phalanx
- Accommodates a Precision<sup>®</sup> Guide cross screw outside the plate which balances the construct to limit plantar gapping
- Plates accommodate Gorilla<sup>®</sup> R3CON™ 2.7, 3.5, and 4.2 mm locking and non-locking polyaxial screws in any plate hole
- Instrumentation compliments the system to address joint preparation and compression at the arthrodesis site




DORSIFLEXION



**PRECISION<sup>®</sup> GUIDE MTP**

- Patent pending guide for trajectory of cross-screw that attaches directly to plate
- Designed to miss all other screws in the construct
- Provides five trajectories of wire paths for variations among patient anatomy
- Places screw outside the plate balancing the construct in all planes
- Allows for 3.0, 3.5, or 4.0 mm cannulated cross screw

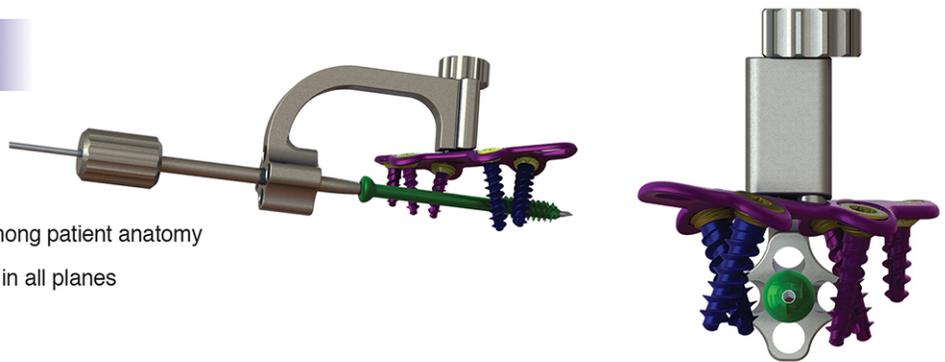

**Gorilla<sup>®</sup> MTP Plates (32 Total Plates)**



Short



Primary



Revision



Graft-Spanning

|                          | Short                               | Primary   | Revision   | Graft-Spanning             |
|--------------------------|-------------------------------------|---|--|----------------------------|
| Total Plates             | 6                                   | 18  | 6  | 2                          |
| Plate Thickness          | 1.3 mm                              | 1.3 mm  | 1.6 mm   | 1.6 mm                     |
| Hole Configurations      | 4 Hole w/ Compression Slot          | 5 Hole w/ Compression Slot<br>6 Hole w/ Compression Slot          | 6 Hole w/ Compression Slot<br>7 Hole w/ Compression Slot<br>9 Hole w/ Compression Slot | 7 Hole w/ Compression Slot |
| Size                     | 0 degree,<br>5 degree,<br>10 degree | 0 degree (S, M, L);<br>5 degree (S, M, L);<br>10 degree (S, M, L) | 5 degree (S, M, L)   | 5 degree                   |
| Left and Right Offerings | Yes                                 | Yes   | Yes  | Yes                        |

# Gorilla® MTP Plating System

## ADDITIONAL INFORMATION




### SPIN GUARD REAMERS

- Available in 4 diameters: 17, 19, 21 and 23 mm
- Comes standard in every PRECISION® Guide MTP Caddy
- Sleeve designed to protect adjacent bone during cup and cone reaming




### PRESERVE® MTP ALLOGRAFT

- Donor harvest site is density matched specific to the MTP indication for strength demands and blood flow requirements
- Aseptically processed without gamma irradiation or hydrogen peroxide to help preserve the native mechanical advantages of human bone and the osteoinductivity of the environment from which the graft is being harvested
- 8 unique sizes available to account for varying anatomies
- Ball and socket created with graft allows precise adjustment in hallux placement by allowing for positioning in all three planes
- Grafts are cut with the exact same reamers as are included in the Gorilla® MTP Plating System ensuring the graft matches the reaming done by surgeon intraoperatively
  - Restores length to the first ray
- Graft is intended to be used in combination with Gorilla® MTP Graft Spanning Plate which plates dorsally and allows for a compression cross screw to be placed through the graft

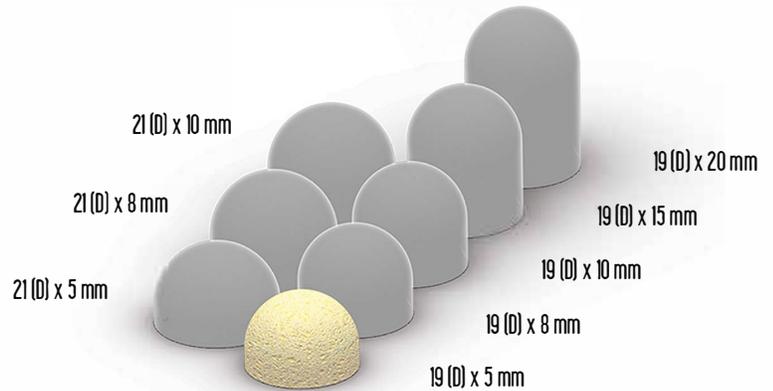

### MTP Trial Sizers

- Allows surgeon to determine ideal graft size by demonstrating correction

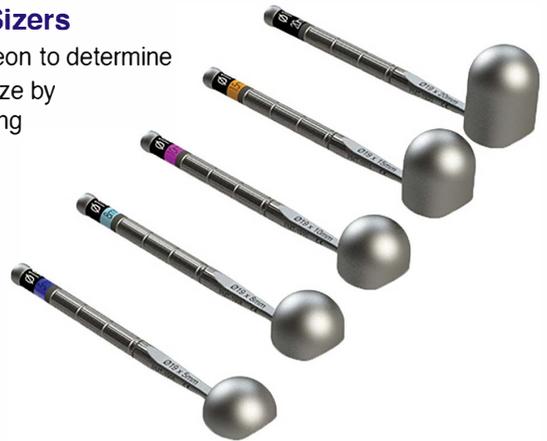

For the contraindications, potential complications and adverse reactions, warnings and precautions associated with this device, please refer to the device specific instructions for use at <http://www.paragon28.com/ifus>

MPSS-01 RevE

™ Trademarks and ® Registered Trademarks of Paragon 28, Inc.

© Copyright 2020 Paragon 28, Inc. All rights reserved.

Patents: [www.paragon28.com/index.php/patents](http://www.paragon28.com/index.php/patents)

Paragon 28, Inc.

14445 Grasslands Dr.  
Englewood, CO 80112  
(855) 786-2828



Paragon 28 Medical Devices Trading Limited

First Floor Block 7 Beckett Way

Park West Business Park

Dublin 12

D12 X884

Ireland

+353 (0) 1588 0350

Exclusively foot & ankle Paragon®