



LISFRANC PLATING SYSTEM


## PRODUCT INFORMATION







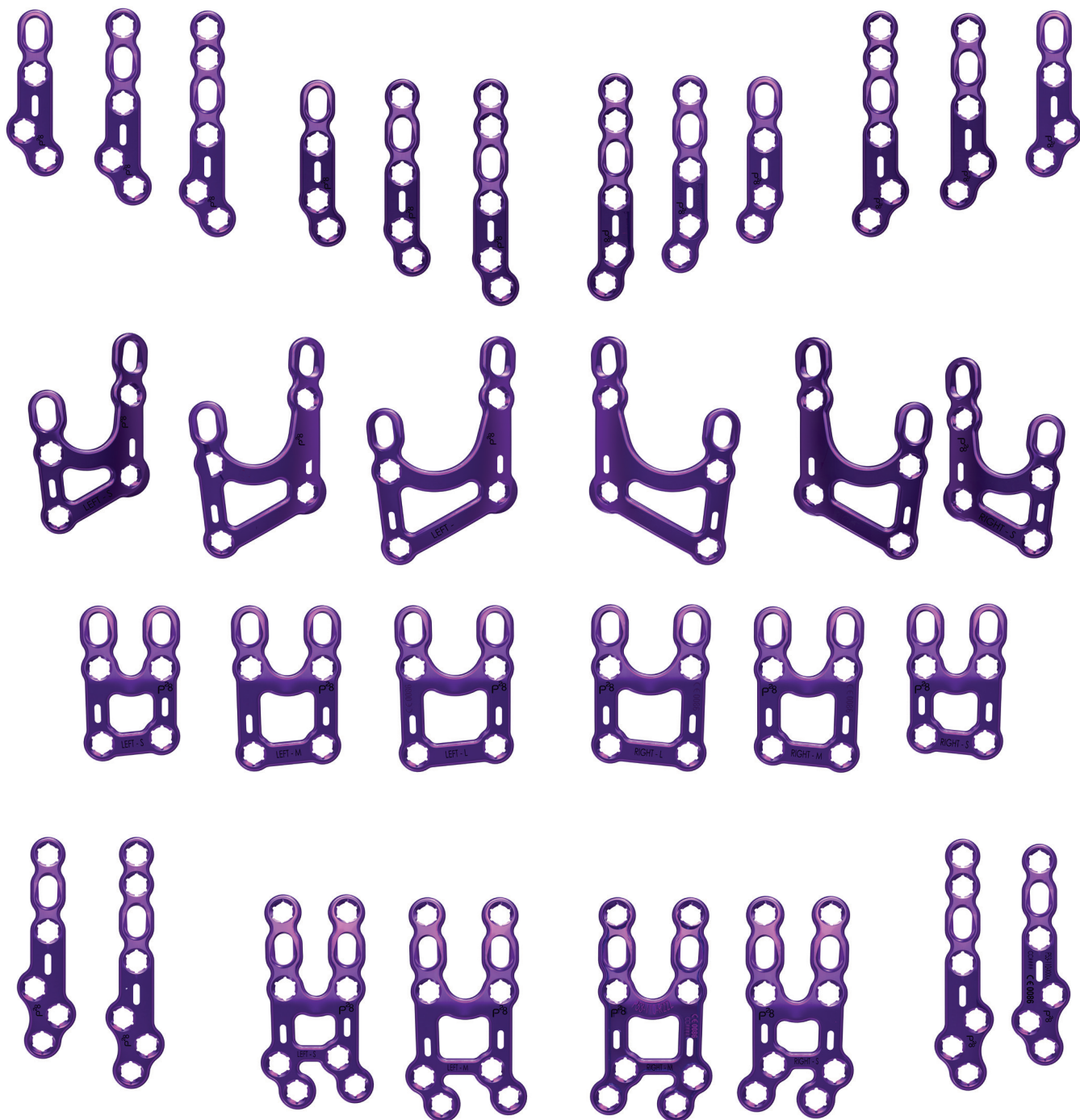


Exclusively foot & ankle <sup>20</sup>  
**Paragon**®

#### FEATURES & BENEFITS

- 6 Unique plating styles
  - ▶ 32 Plating options
- 1.4 mm low profile plate design
- Anatomically contoured for TMT joints




## LISFRAC/TMT PLATING OPTIONS

	Dual Ray			Single Ray		
	1 <sup>st</sup> & 2 <sup>nd</sup>	2 <sup>nd</sup> & 3 <sup>rd</sup>		Slanted Straight	Slanted T	Clover
						
<b>Total Plates</b>	6	4	6	6	6	6
<b>Hole Configurations</b>	4 Hole w/ 2 Compression Slots Small, Medium, Large	8 Hole w/ 2 Compression Slots Small, Medium	4 Hole w/ 2 Compression Slots Small, Medium, Large	3 Hole w/ Compression Slot 4 Hole w/ Compression Slot 5 Hole w/ Compression Slot	3 Hole w/ Compression Slot 4 Hole w/ Compression Slot 5 Hole w/ Compression Slot	5 Hole w/ Compression Slot 6 Hole w/ Compression Slot

## PLATE FEATURES

- 1.4 mm plate thickness helps to reduce soft tissue irritation
- Made from titanium alloy Ti 6Al-4V ELI allowing for increased plate strength
- Anatomic contour can help reduce intraoperative bending which may save time in OR
- All plates contain compression slot which sits distally allowing for better screw purchase in diaphyseal bone
- Hook slot for use with Compression/Distractor Device allowing in-line compression with the plate
- All holes accommodate variable angle locking screw technology (up to 15 degrees off-axis in any direction) to allow for targeting of better quality bone and avoidance of the joint space

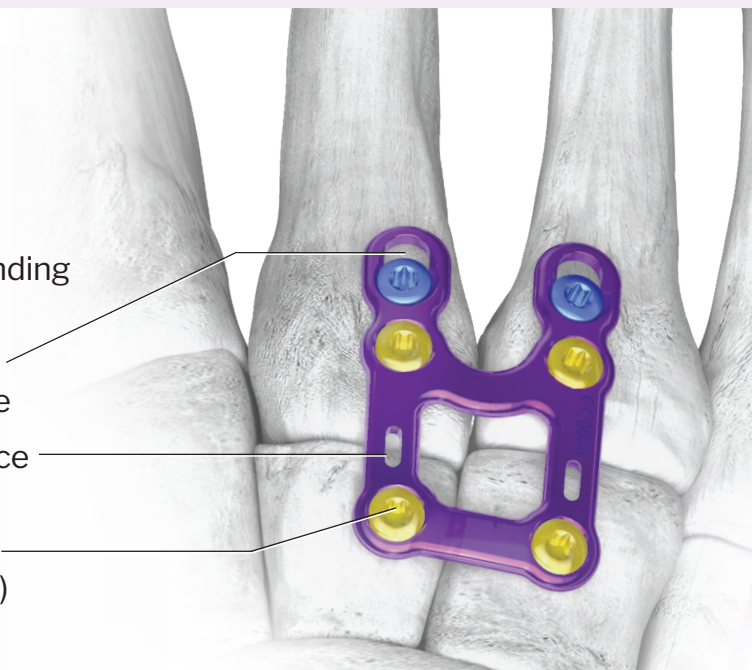



ADDITIONAL PARAGON 28<sup>®</sup> SOLUTIONS



T-Plate  
8 Hole  
Compression T2



L-Plate  
6 Hole  
Compression L2  
Left




Oblique T-Plate  
8 Hole  
Compression T2  
Left



Oblique L-Plate  
7 Hole  
Compression  
Left



Y-Plate  
7 Hole  
Compression Y3  
Left

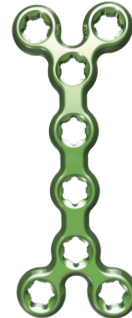


Y-Plate  
8 Hole  
Dual Y 3-2-3



2.5 mm Mini-Monster<sup>®</sup>  
Solid Screw



3.0 mm Mini-Monster<sup>®</sup>  
Solid Screw

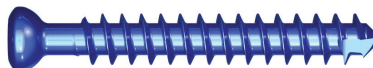

4.0 mm Mini-Monster<sup>®</sup>  
Solid Screw



LSFR-01 Rev D

2021-07-29

™Trademarks and ®Registered Marks of Paragon 28<sup>®</sup>, Inc.

© Copyright 2021 Paragon 28<sup>®</sup>, Inc. All rights reserved.

Patents: [www.paragon28.com/patents](http://www.paragon28.com/patents)

Paragon 28, Inc. 

14445 Grasslands Dr.

Englewood, CO 80112 USA

(855) 786-2828

Paragon 28 Medical Devices Trading Limited

First Floor Block 7 Beckett Way

Park West Business Park

Dublin 12 D12 X884

Ireland

+353 (0) 1588 0350



[www.Paragon28.com](http://www.Paragon28.com)