


R3LEASE™

STABILIZATION SYSTEM



R3LEASE™ Stabilization System

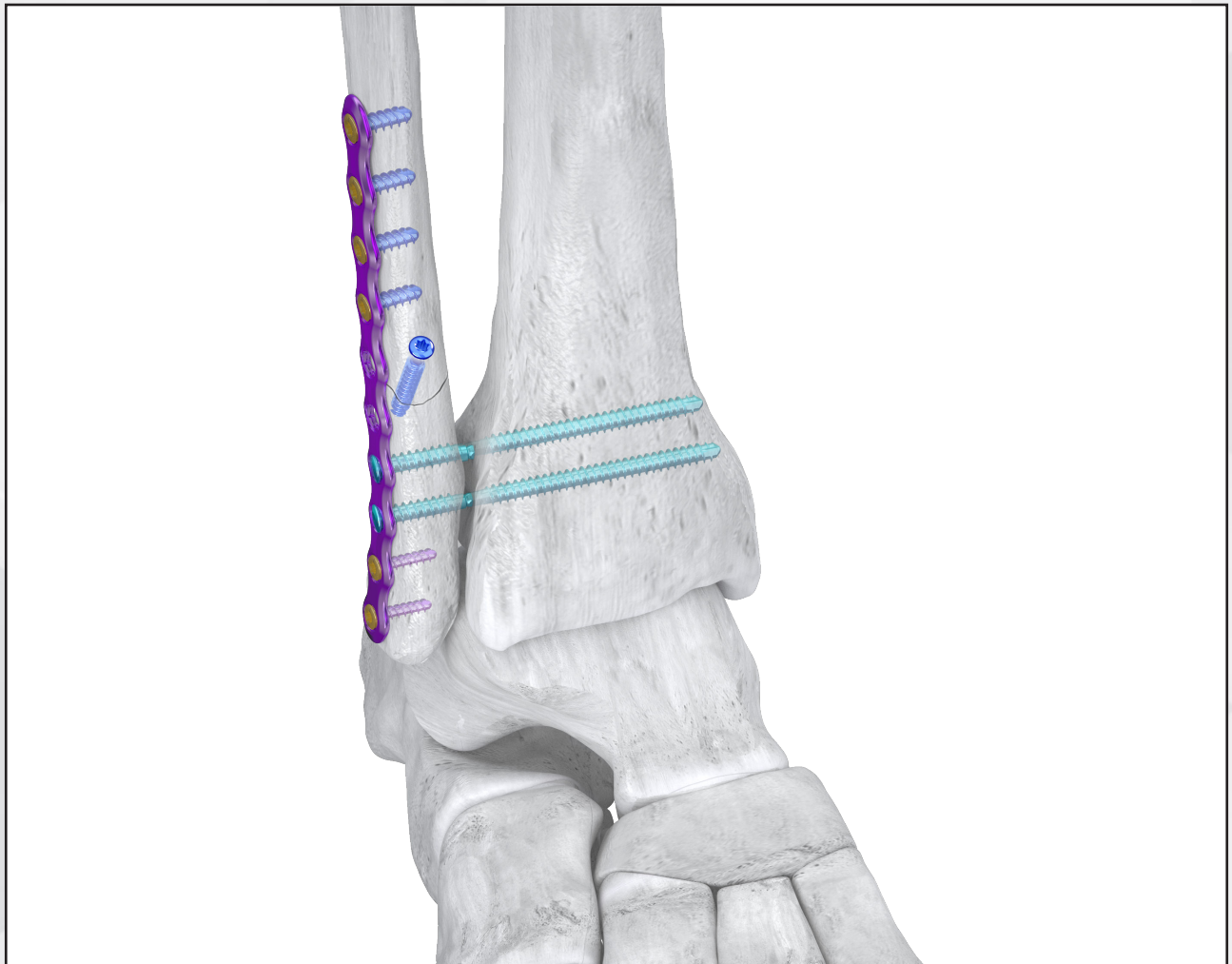


**Shown with the Gorilla® R3CON Plating System*

PRODUCT INFORMATION

Designed to Provide a Simple Solution for Syndesmosis Repair


- ▶ Once a patient begins weightbearing, traditional syndesmotic screws may fracture
- ▶ The location of screw fracture of a traditional syndesmotic screw can vary, leading to difficult removal situations and/or potential pain at the screw fracture site
- ▶ The R3LEASE™ Stabilization System was designed so that if a screw fracture occurs, the screw breaks cleanly at the notch point in the clear space

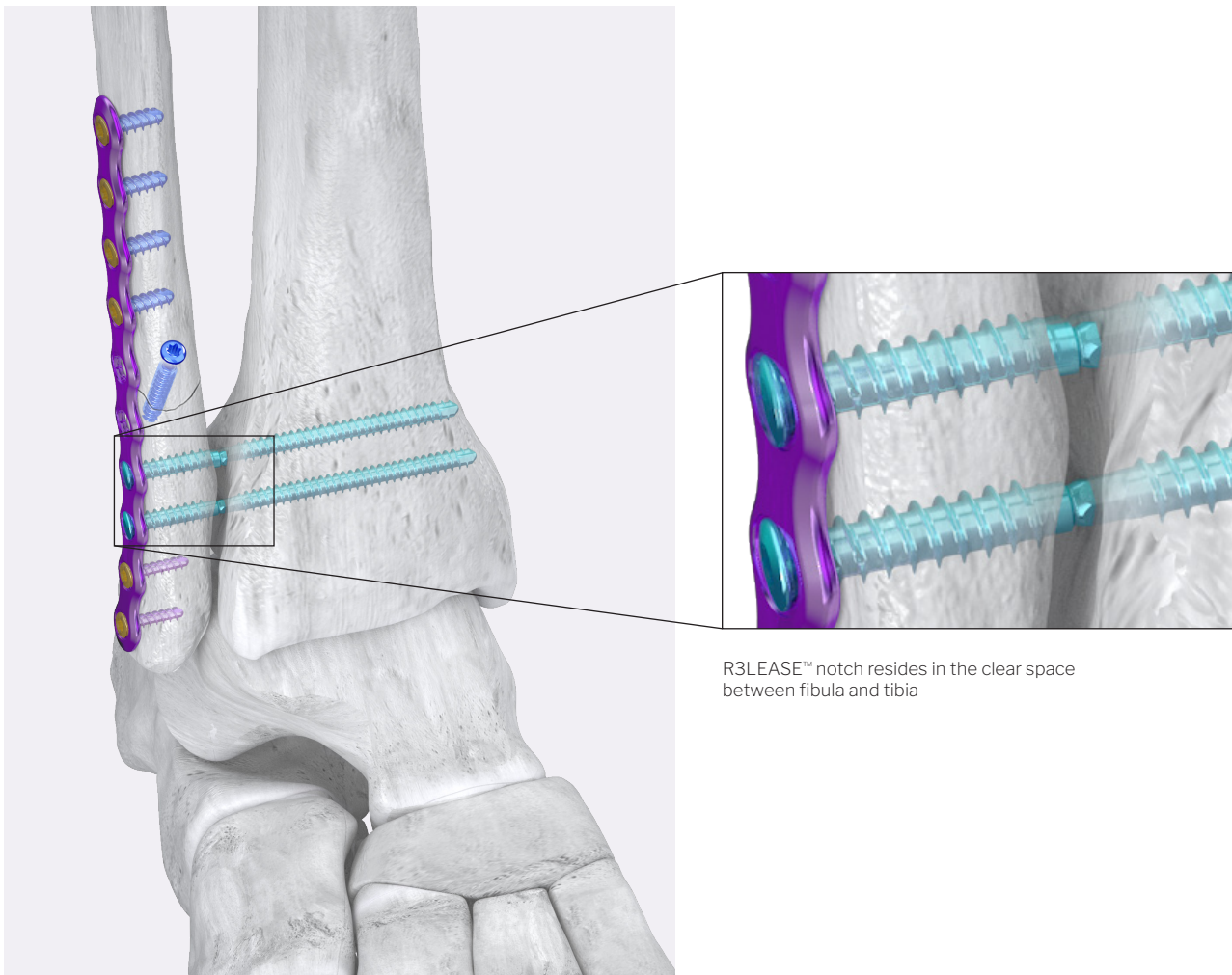

R3LEASE™ notch resides in the clear space between fibula and tibia


R3LEASE™ SCREW SIZES

Ø3.9 mm R3LEASE™ Screw 14 mm Notch Length

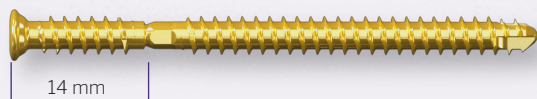

Screw Lengths:
40 mm–45 mm (in 5 mm increments)
48 mm–64 mm (in 2 mm increments)



Ø3.9 mm R3LEASE™ Screw 17 mm Notch Length

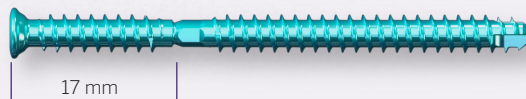

Screw Lengths:
40 mm–45 mm (in 5 mm increments)
48 mm–64 mm (in 2 mm increments)

SCREW FEATURES

- ▶ Multiple Notch Lengths: 14 mm and 17 mm
 - Addresses differences in patient anatomy and accommodates use with or without a plate
 - 2 mm increments from 48–64 mm intended for tricortical or quadricortical fixation
- ▶ Strength in Design
 - Static bending strength and rigidity of the R3LEASE™ Screw exceeded that of a comparable 3.5 mm solid screw as well as a 4.0 mm cannulated screw made of the same material.¹
- ▶ If removal is desired, multiple removal features have been incorporated:

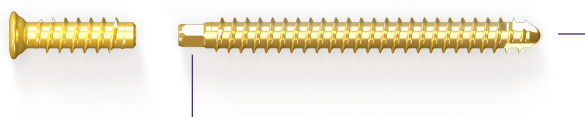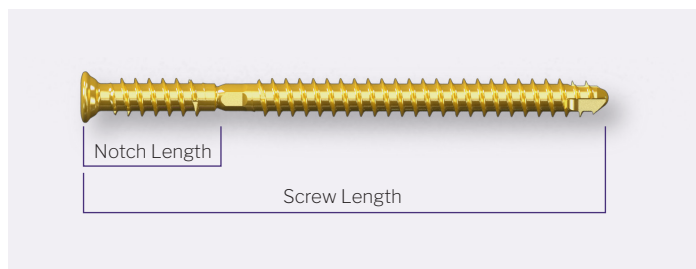

- ▶ Lateral removal feature
 - Hex feature allows for lateral removal of tibial fragment through the fibula
- ▶ Medial removal feature
 - Medial removal feature on screw tip allows for medial removal using a specialized driver

Washer



Washers available for R3LEASE™ screws

R3LEASE™

STABILIZATION SYSTEM




R3LEASE™ Stabilization System

Reference:

1. Internal data on file, TR-17102601

GBSS-01 RevE
2023-10-12

Paragon 28°, Inc. 
14445 Grasslands Dr.
Englewood, CO 80112 USA
(855) 786-2828

Paragon 28° Medical Devices Trading Limited
First Floor Block 7 Beckett Way
Park West Business Park
Dublin 12
D12 X884
Ireland
+353 (0) 1588 0350

Exclusively foot & ankle ²⁸
Paragon®

www.Paragon28.com

™Trademarks and *Registered Marks of Paragon 28°, Inc.
© Copyright 2023 Paragon 28°, Inc. All rights reserved.
Patents: www.paragon28.com/patents

For the contraindications, potential complications and adverse reactions, warnings and precautions associated with this device, please refer to the device specific instructions for use at <http://www.paragon28.com/ifus>