





Gorilla® Lateral Column Plating System

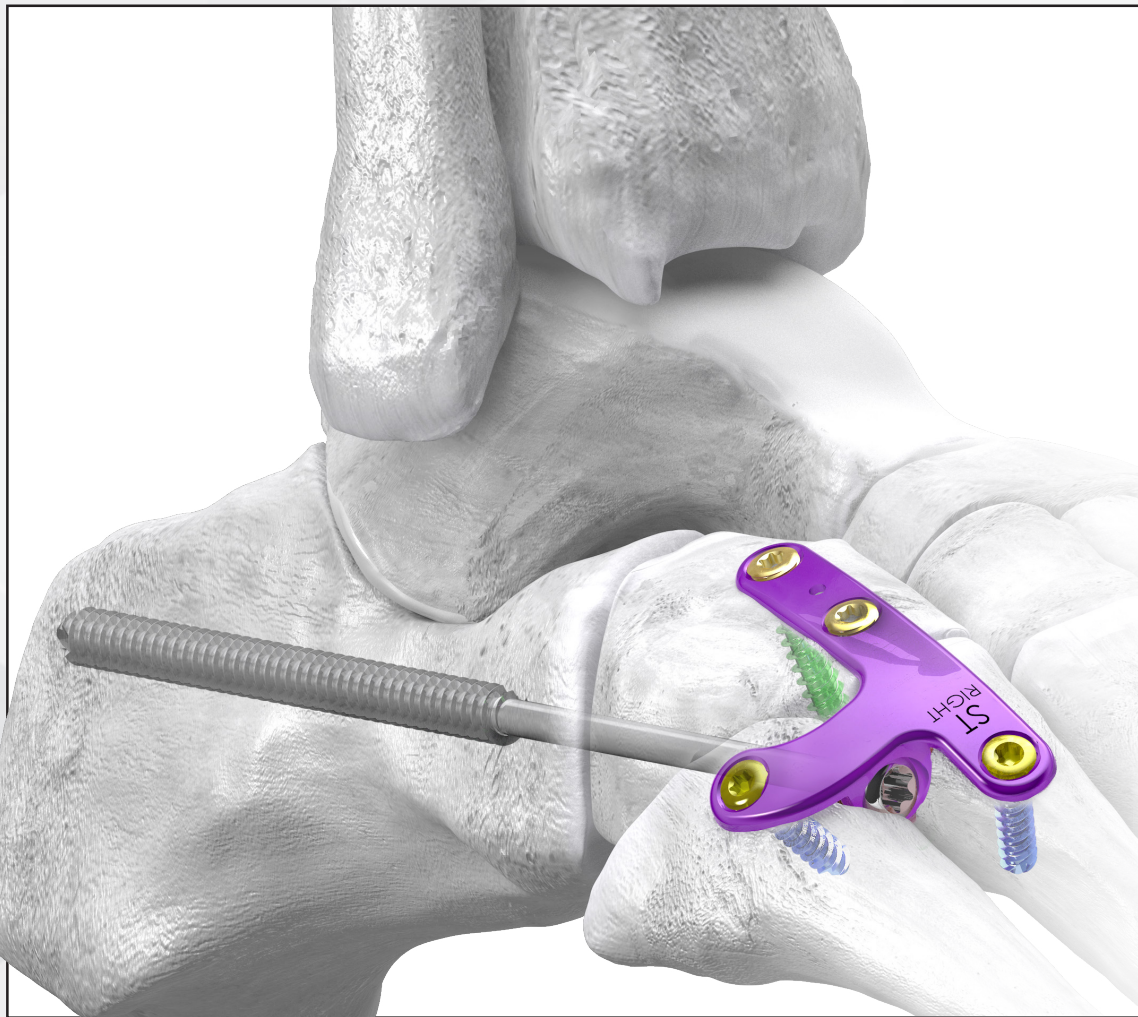

GORILLA® LATERAL COLUMN PLATING SYSTEM

Features and Benefits:

- ▶ In the case where the cuboid is the cause of plantar ulcerations, anatomic reduction of the cuboid may help alleviate plantar pressure and aid in the healing of plantar ulcerations
- ▶ Plate and screw construct was designed to maintain anatomic alignment of the lateral column and prevent plantar subluxation of the cuboid
- ▶ The Gorilla® style Lateral Column Plate is designed to accept a Type II Anodized 5.5 mm beaming plate screw to help stabilize and compress the lateral column

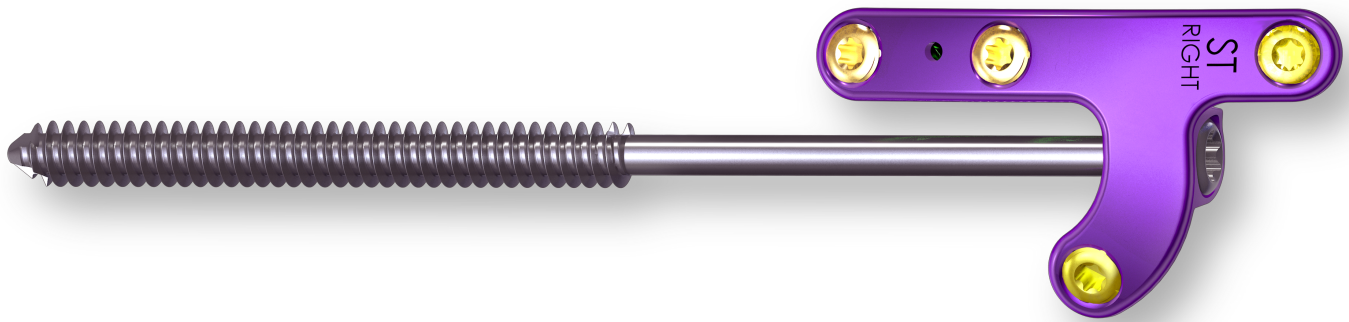

Lateral Column Plate

- ▶ Offered in standard and large sizes
- ▶ 1.5 mm plate thickness
- ▶ Anatomically contoured

Standard plate screw holes accept 2.7 mm, 3.5 mm and 4.2 mm locking and non-locking variable angle plate screws


Proximal screws on the cuboid help maintain anatomical alignment of the lateral column



Recessed screw hole accepts a Type II Anodized 5.5 mm beaming plate screw

Type II Anodized 5.5 mm Beaming Plate Screw

- ▶ Offered in fully threaded and partially threaded options
- ▶ 5.5 mm screw interfaces with plate to help keep the cuboid from subluxing plantarly
- ▶ Offered in 80 mm–120 mm lengths by 5 mm increments

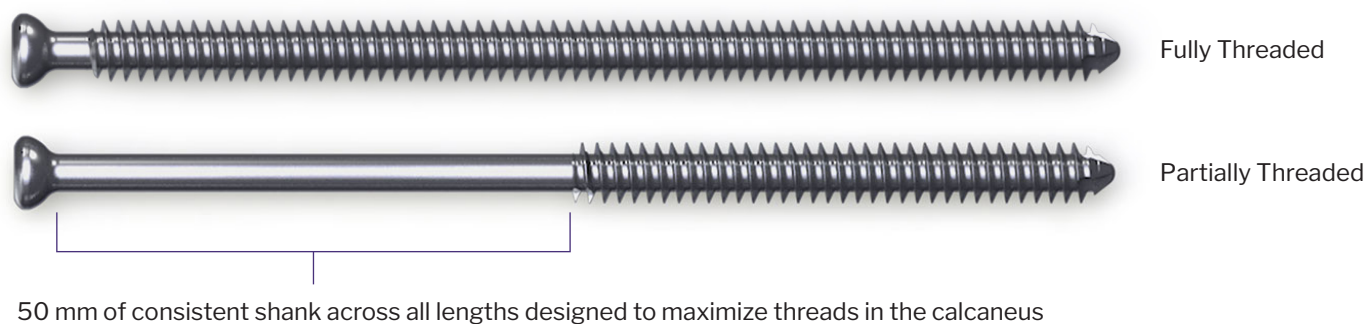

SYSTEM SPECIFIC INSTRUMENTATION

Templating System

- ▶ Designed to establish correct reaming location between the 4th and 5th metatarsals to accommodate space for the 5.5 mm plate screw
- ▶ Joint indicating tab aligns template at the metatarsal-cuboid joint
- ▶ Wire holes through the template help facilitate easy transfer from template to implant



Joint Preparation Paddle

- ▶ Helps center the template in the inter-digital space so equal reaming of the 4th and 5th metatarsals can be achieved




Reamer

- ▶ Laser markings on reamer indicate appropriate reaming depth to accommodate recessed plate hole



Gorilla® Lateral Column Plating System

GISS-01 Rev E
2023-09-05

Paragon 28®, Inc. 
14445 Grasslands Dr.
Englewood, CO 80112 USA
(855) 786-2828

Paragon 28® Medical Devices Trading Limited
First Floor Block 7 Beckett Way
Park West Business Park
Dublin 12
D12 X884
Ireland
+353 (0) 1588 0350

Exclusively foot & ankle ²⁸
Paragon®

www.Paragon28.com

™Trademarks and ®Registered Marks of Paragon 28®, Inc.
© Copyright 2023 Paragon 28®, Inc. All rights reserved.
Patents: www.paragon28.com/patents

For the contraindications, potential complications and adverse reactions, warnings and precautions associated with this device, please refer to the device specific instructions for use at <http://www.paragon28.com/ifus>