




Gorilla® Lisfranc Plating System



PLATING SYSTEM FEATURES

- ▶ 6 Unique plating styles
 - 32 Plating options
- ▶ 1.4 mm low profile plate design
- ▶ Anatomically contoured for TMT joints

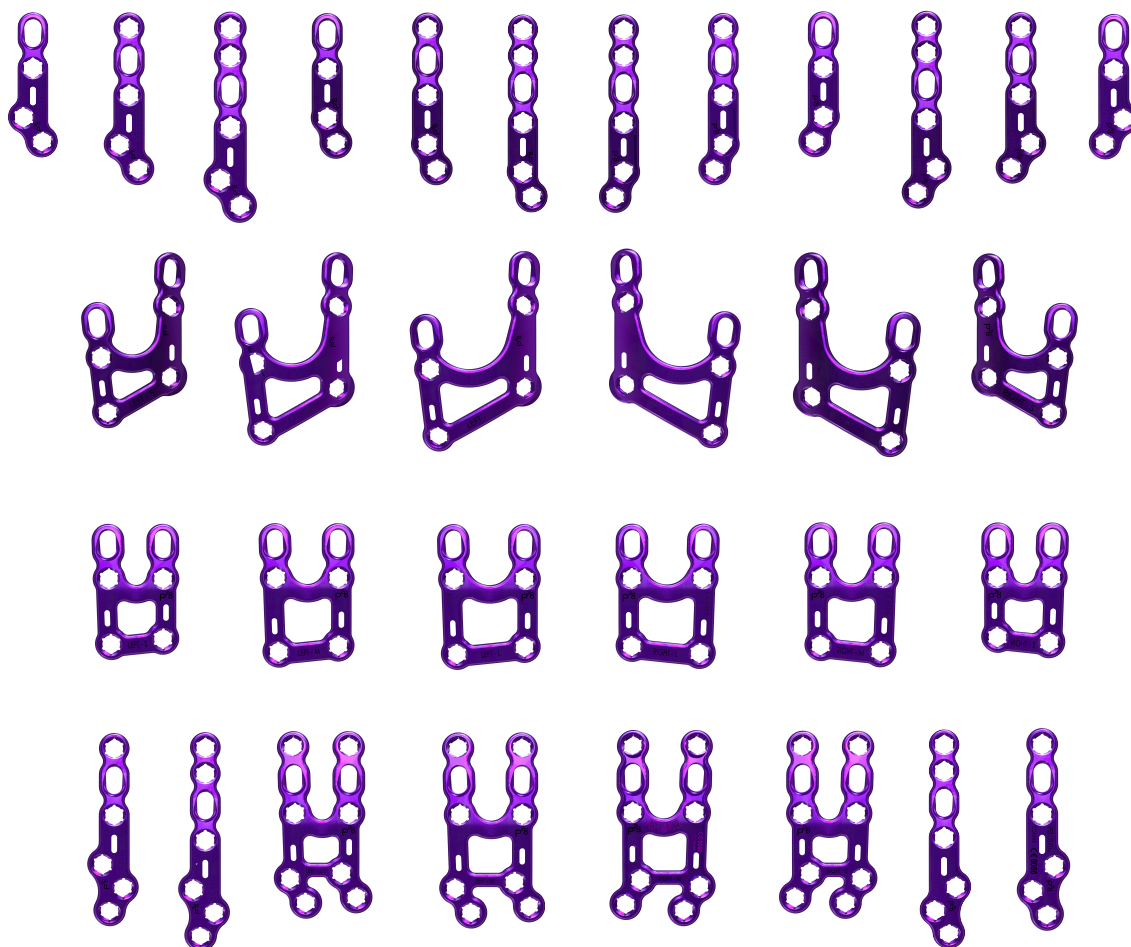








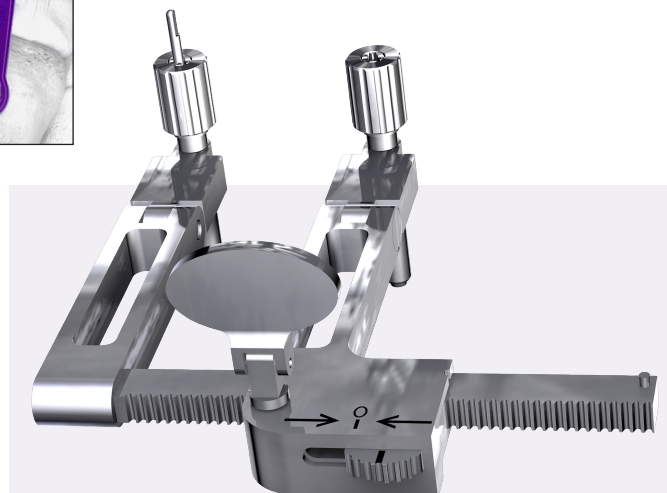
	DUAL RAY			SINGLE RAY		
	1 st & 2 nd	2 nd & 3 rd		Slanted Straight	Slanted T	Clover
						
Total Plates	6	4	6	6	6	6
Hole Configurations	4-Hole w/ 2 Compression Slots Small, Medium, Large	8-Hole w/ 2 Compression Slots Small, Medium	4-Hole w/ 2 Compression Slots Small, Medium, Large	3-Hole w/ Compression Slot 4-Hole w/ Compression Slot 5-Hole w/ Compression Slot	3-Hole w/ Compression Slot 4-Hole w/ Compression Slot 5-Hole w/ Compression Slot	5-Hole w/ Compression Slot 6-Hole w/ Compression Slot

PLATE FEATURES

- ▶ 1.4 mm plate thickness helps to reduce soft tissue irritation
- ▶ Made from titanium alloy Ti 6Al-4V ELI allowing for increased plate strength
- ▶ Anatomic contour can help reduce intraoperative bending which may save time in OR
- ▶ All plates contain compression slot which sits distally allowing for better screw purchase in diaphyseal bone
- ▶ Hook slot for use with Compression/ Distraction Device allowing in-line compression with the plate
- ▶ All holes accommodate variable angle locking screw technology (up to 15 degrees off-axis in any direction) to allow for targeting of better quality bone and avoidance of the joint space

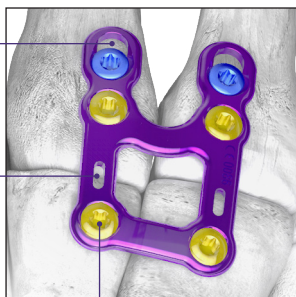

Compression/Distraction Device

ADDITIONAL PARAGON 28® SOLUTIONS

Baby Gorilla® Mini Plating System



T-Plate
8-Hole
Compression T2



L-Plate
6-Hole
Compression L2
Left



Oblique T-Plate
8-Hole
Compression T2
Left



Oblique L-Plate
7-Hole
Compression
Left




Y-Plate
7-Hole
Compression Y3
Left




Y-Plate
8-Hole
Dual Y 3-2-3


Mini-Monster® Screw System

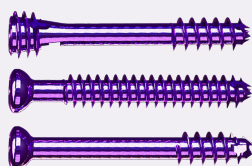

Ø3.0 mm Mini-Monster®
Solid Screw

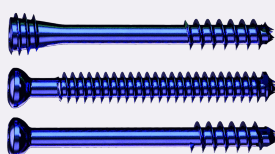

Ø3.5 mm Mini-Monster®
Solid Screw

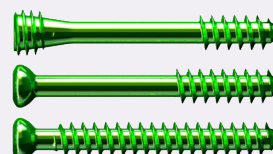


Ø4.0 mm Mini-Monster®
Solid Screw



Gorilla® Lisfranc Plating System

LSFR-01 RevF
2023-10-27

Paragon 28®, Inc. 
14445 Grasslands Dr.
Englewood, CO 80112 USA
(855) 786-2828

Paragon 28® Medical Devices Trading Limited
First Floor Block 7 Beckett Way
Park West Business Park
Dublin 12
D12 X884
Ireland
+353 (0) 1588 0350

™Trademarks and ®Registered Marks of Paragon 28®, Inc.
© Copyright 2023 Paragon 28®, Inc. All rights reserved.
Patents: www.paragon28.com/patents

Exclusively foot & ankle ²⁸
Paragon®

www.Paragon28.com

For the contraindications, potential complications and adverse reactions, warnings and precautions associated with this device, please refer to the device specific instructions for use at <http://www.paragon28.com/ifus>