



PARATROOPER[®]

PLANTAR PLATE REPAIR SYSTEM

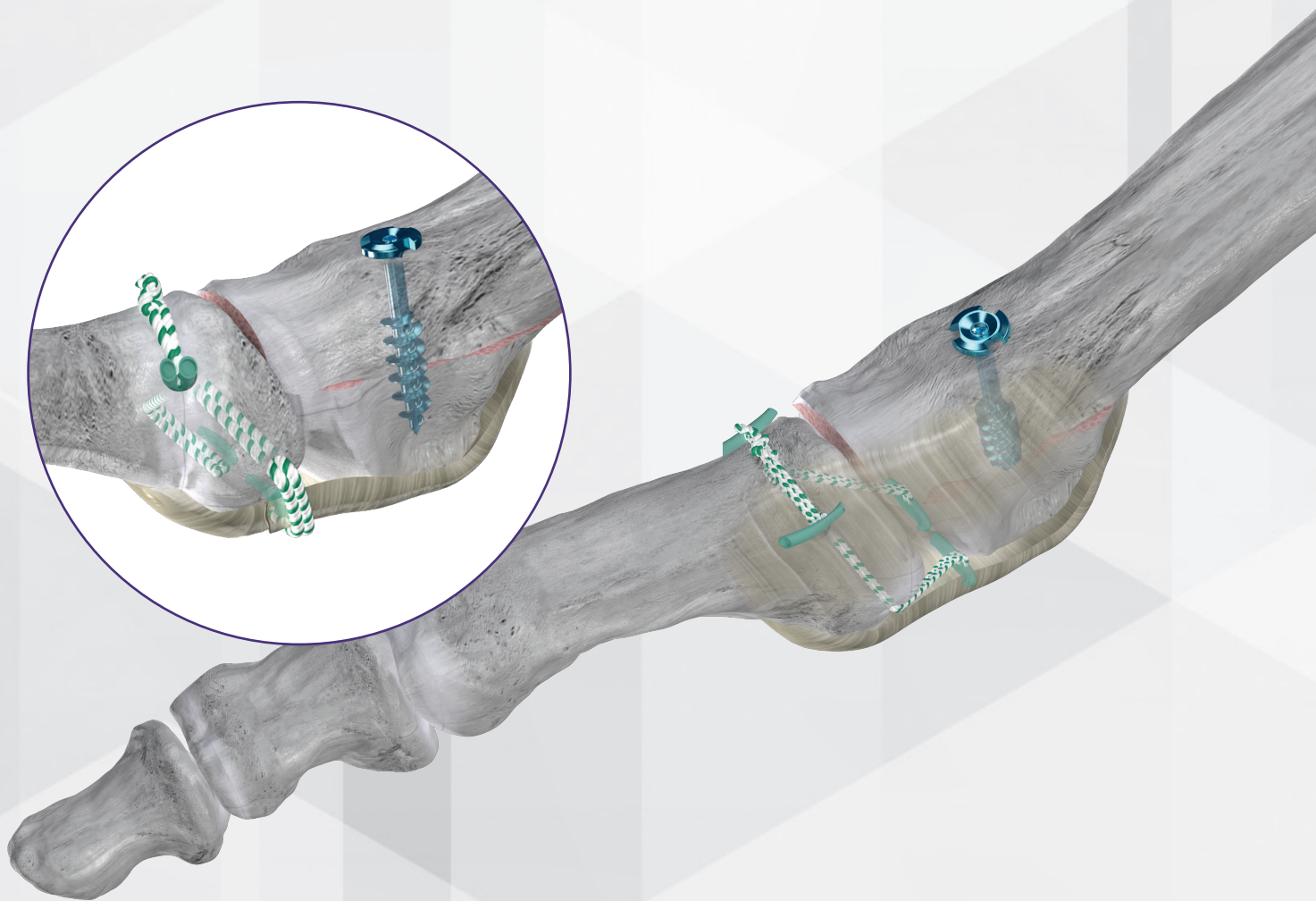

ALL-SUTURE IMPLANT

- ▶ Designed for optimal insertion into tissue or bone
- ▶ Slight tension after insertion will ball the suture “sock”, locking the implant in place
 - Can increase or release tension until knots are tied
- ▶ The implant features two connected all-suture anchors
 - Provides for independent tensioning with single or dual implant construct
 - Differential tightening assists with varus and valgus correction
- ▶ Implant fits through a small 1.4 mm bone tunnel
 - *Comes attached to plantar plate specific needle to provide consistent and reliable delivery*

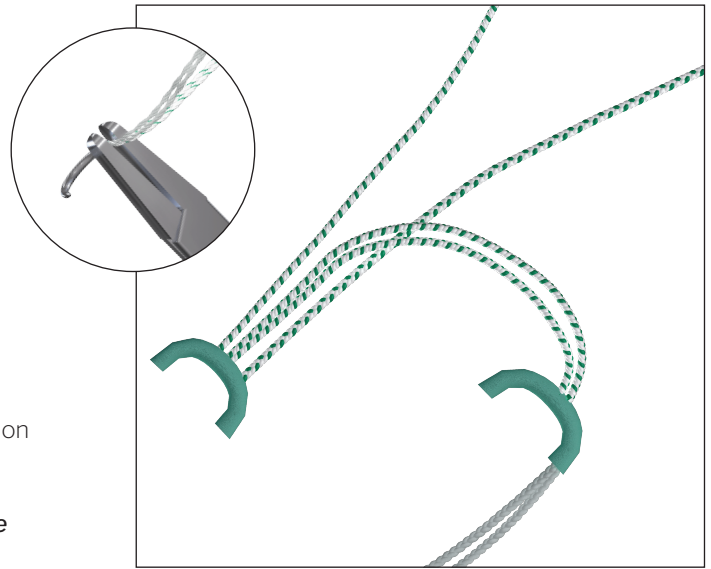


ALL IN ONE – PLANTAR PLATE REPAIR

One Implant For Every Approach


- ▶ System accommodates either a plantar or dorsal approach
 - Both utilize the same sequentially tensioned implant system
- ▶ Weil Kit with Monster® BITE screws for dorsal approach or concurrent Weil osteotomies
 - *Monster® BITE screws feature a unique torque limit designed to break the pin away when the screw head contacts the near cortex of bone*




DORSAL APPROACH
Plate to Bone

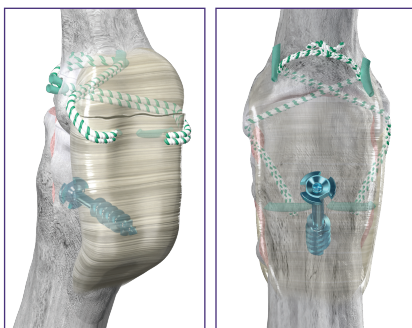


PLANTAR APPROACH
Plate to Bone

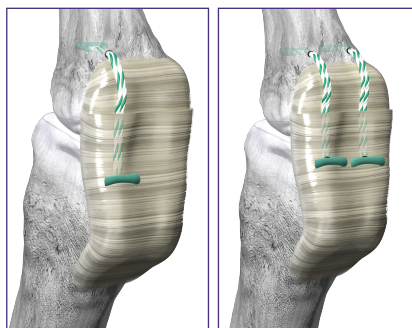

PLANTAR APPROACH
Plate to Plate

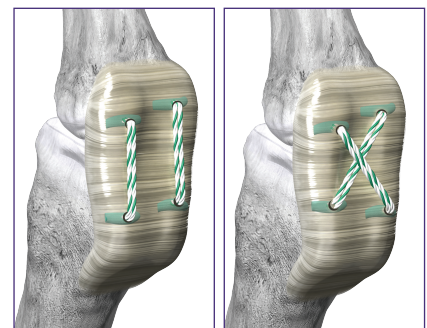

PARATROOPER™ IMPLANT

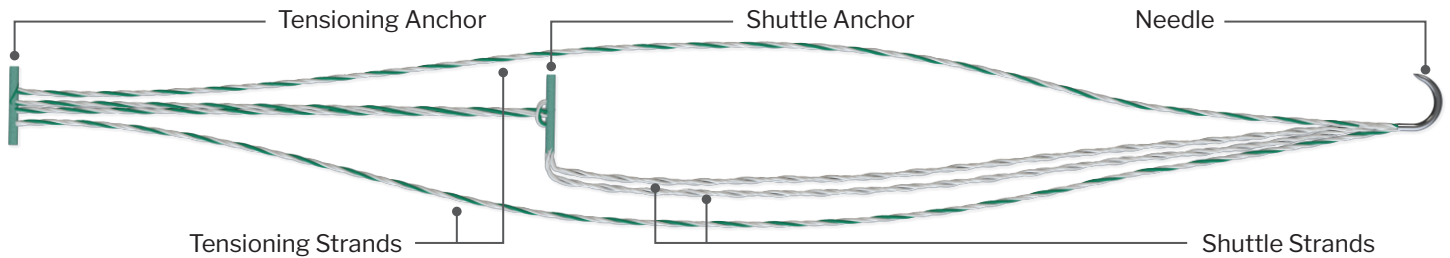

SYSTEM SPECIFIC INSTRUMENTATION

Set instrumentation allows for easy bone preparation and access to the plantar plate using either the dorsal or plantar approach.

Welpi Retractor

Single retractor designed for the plantar approach, eliminates the need for multiple instruments, retracts fat pad and flexors simultaneously. Hook Attachment grabs flexor tendons and attaches to Welpi Retractor.




Welpi Hook Attachments



HY Mini Pin Distractor


Small distractor for exposure to the plantar plate from the dorsal approach.

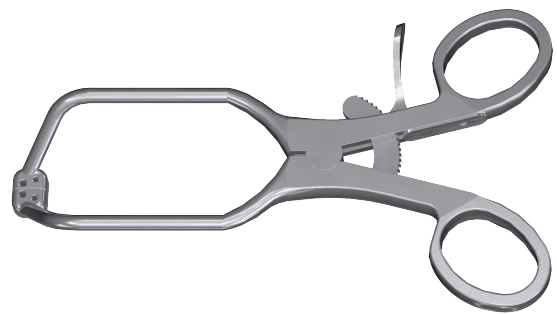


McGlamry Elevator



Weil Osteotomy Clamp


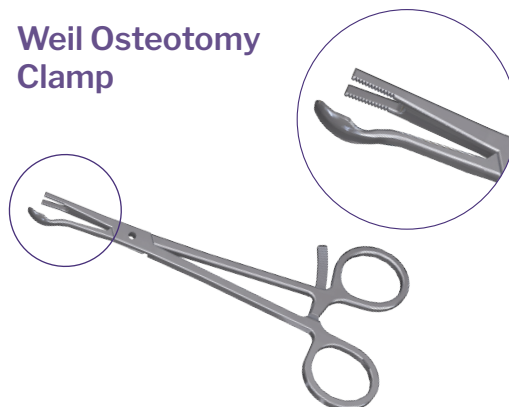


Needle Driver

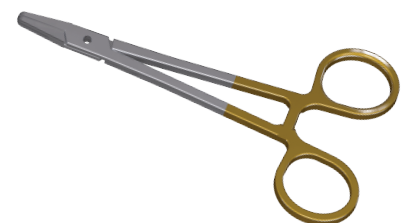

IMPLANT ONLY STERILE KIT

Paratrooper™ Implant

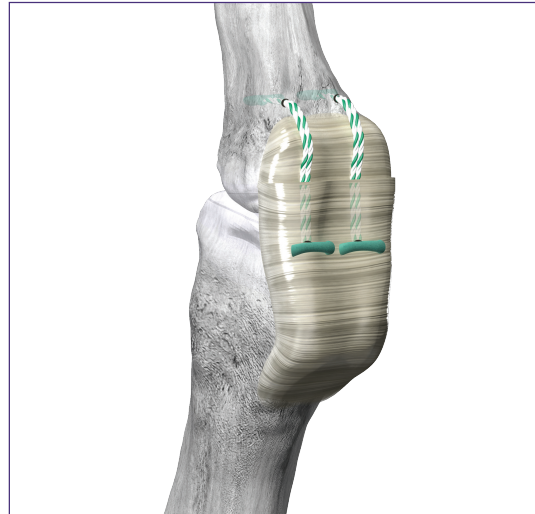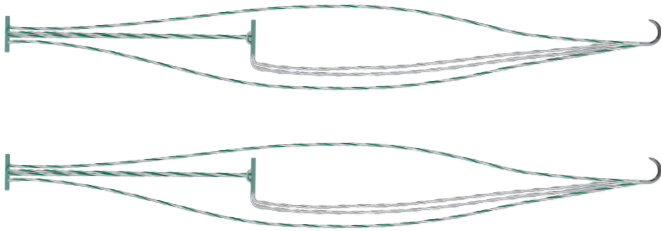


WEIL KIT

Monster® BITE Driver




Monster® BITE Screw Ø2.0 x 11 mm and Ø2.0 x 13 mm

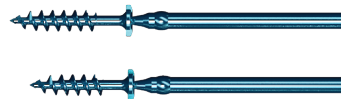


1.4 mm X 15 cm Kirschner Wire



Suture Retriever

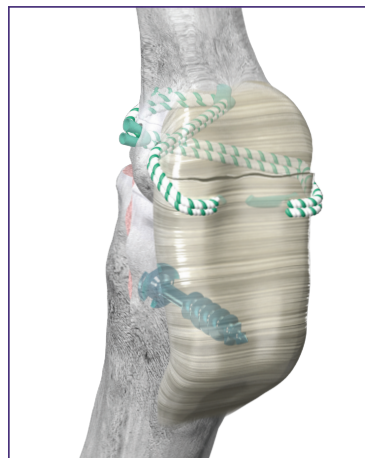

PT-BRO-01 Rev. A
2022-02-28

Paragon 28®, Inc. 
14445 Grasslands Dr.
Englewood, CO 80112 USA
(855) 786-2828

Paragon 28® Medical Devices Trading Limited
First Floor Block 7 Beckett Way
Park West Business Park
Dublin 12
D12 X884
Ireland
+353 (0) 1588 0350

Exclusively foot & ankle ²⁸
Paragon®

www.Paragon28.com

™Trademarks and *Registered Marks of Paragon 28®, Inc.
© Copyright 2022 Paragon 28®, Inc. All rights reserved.
Patents: www.paragon28.com/patents

For the contraindications, potential complications and adverse reactions, warnings and precautions associated with this device, please refer to the device specific instructions for use at <http://www.paragon28.com/ifus>