




TITAN 3-D™ Wedge System

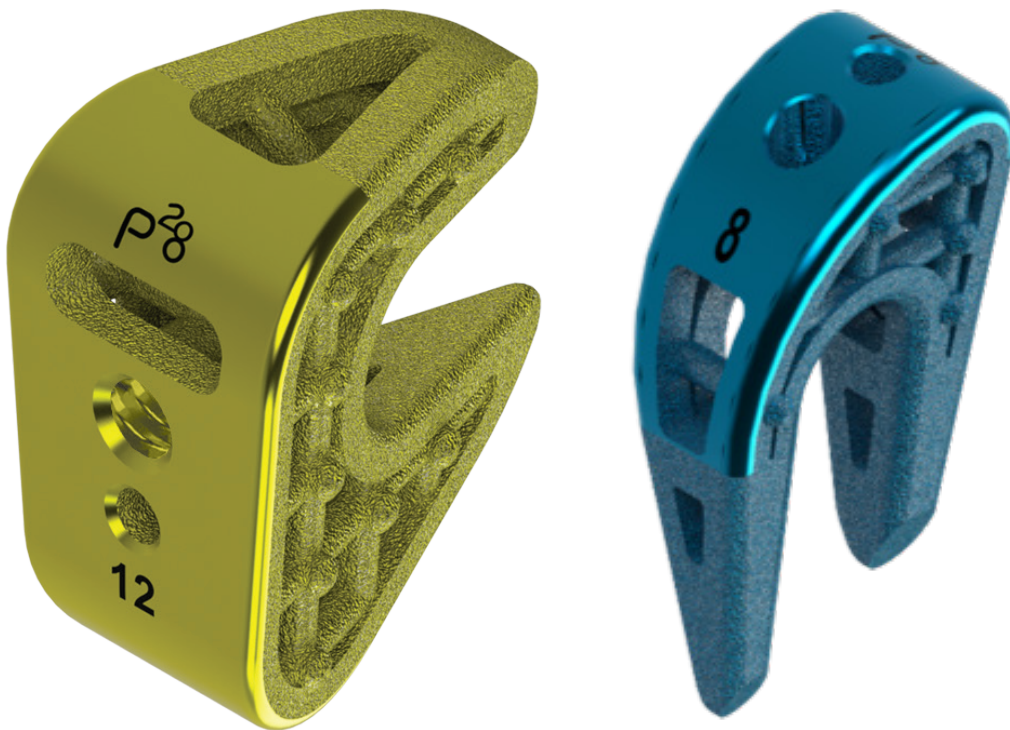

TITAN 3-D™ EVANS WEDGE

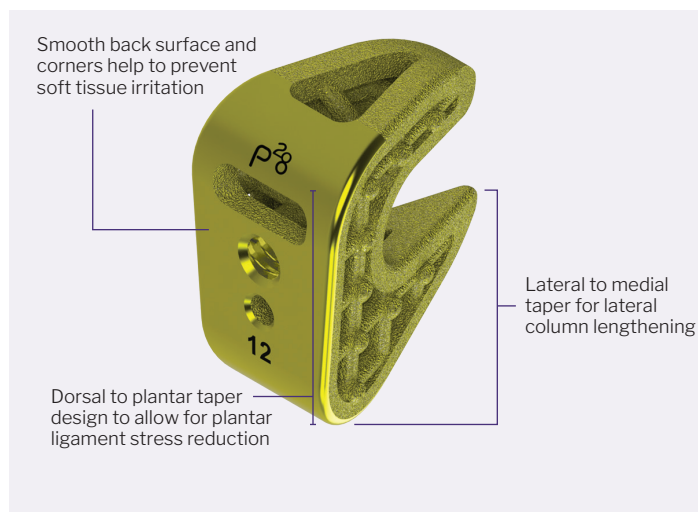


TITAN 3-D™ COTTON WEDGE

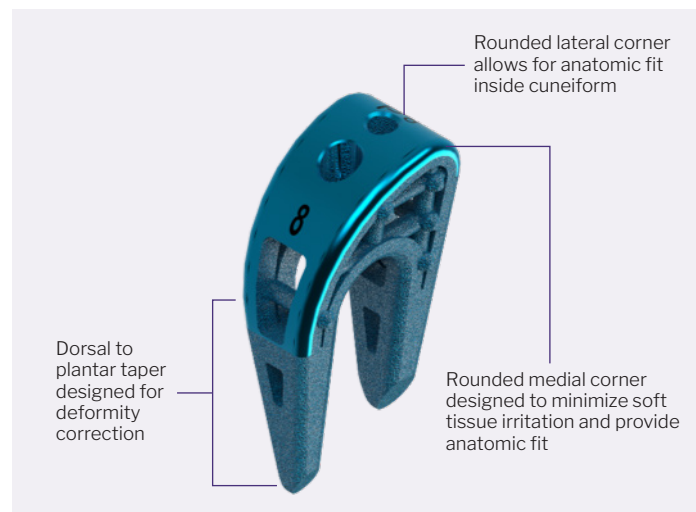



FEATURES PRESENT IN BOTH IMPLANT FAMILIES

- ▶ Manufactured from medical grade titanium alloy (Ti-6Al-4V)
- ▶ Tapered nose helps to aid in implant insertion
- ▶ Open geometry with three-dimensional scaffold allows for blood entry, bone through growth and the incorporation of biologic products, if used
- ▶ Central opening allows for passage of screw across osteotomy to help increase stability of the construct
- ▶ Spikes on both sides of implant designed to help prevent expulsion from osteotomy site
- ▶ Dedicated procedure specific trials to allow surgeon to determine ideal wedge size by demonstrating correction

TITAN 3-D™ WEDGE SYSTEM—12 UNIQUE TITANIUM WEDGES

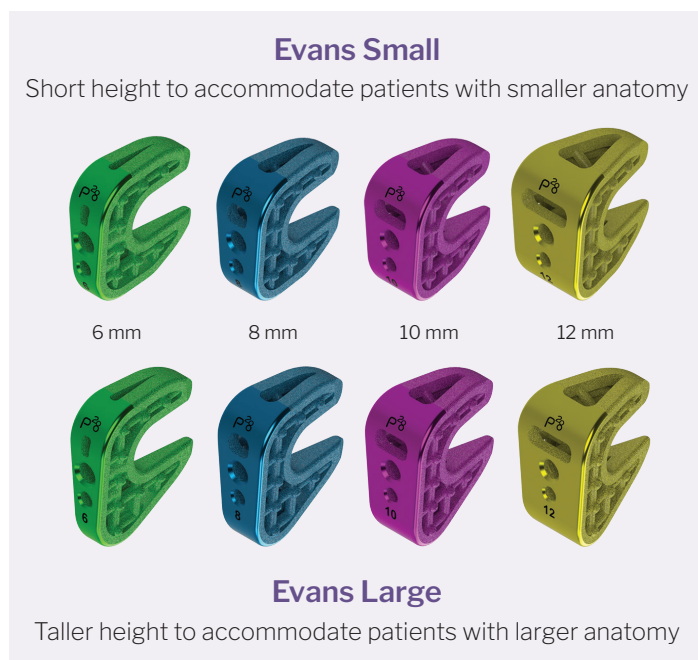


PRODUCT SPECIFIC INSERTER (EVANS AND COTTON INSERTERS)

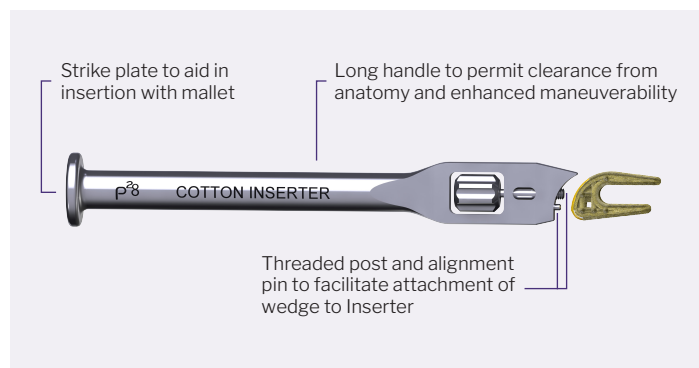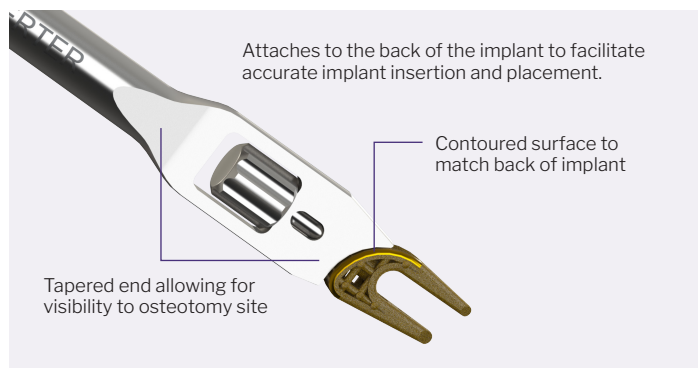


EVANS AND COTTON GUIDES

Precision Guide System

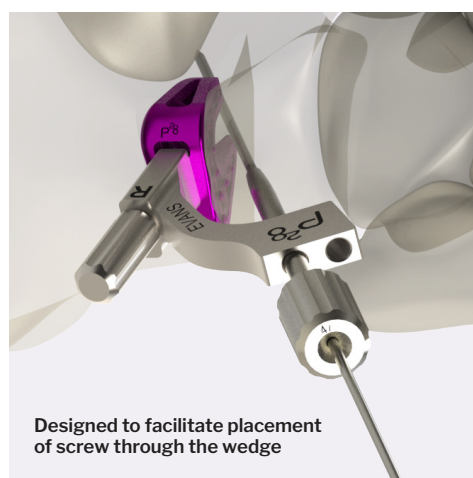


Resection Guides

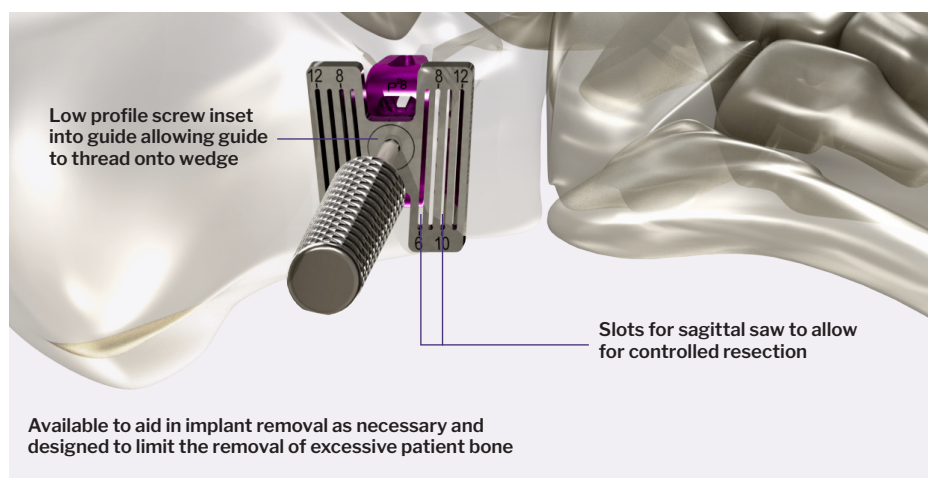

ADDITIONAL FIXATION OPTIONS


(may be inserted through wedge)

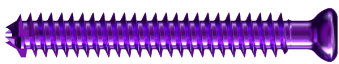Ø3.5 mm and Ø4.0 mm Mini-Monster® Fully Threaded Screws	
Ø3.5 mm Fully Threaded	Ø4.0 mm Fully Threaded
	
10 mm–50 mm (2.0 mm increments)	10 mm–50 mm (2.0 mm increments)
	55 mm–60 mm (5.0 mm increments)



TITAN 3-D™ Wedge System

T3WS-01 RevD
2023-10-12

Paragon 28®, Inc. 
14445 Grasslands Dr.
Englewood, CO 80112 USA
(855) 786-2828

Paragon 28® Medical Devices Trading Limited
First Floor Block 7 Beckett Way
Park West Business Park
Dublin 12
D12 X884
Ireland
+353 (0) 1588 0350

™Trademarks and ®Registered Marks of Paragon 28®, Inc.
© Copyright 2023 Paragon 28®, Inc. All rights reserved.
Patents: www.paragon28.com/patents

Exclusively foot & ankle ²⁸
Paragon®

www.Paragon28.com

For the contraindications, potential complications and adverse reactions, warnings and precautions associated with this device, please refer to the device specific instructions for use at <http://www.paragon28.com/ifus>