




Bun-Yo-Matic™ Lapidus Clamp

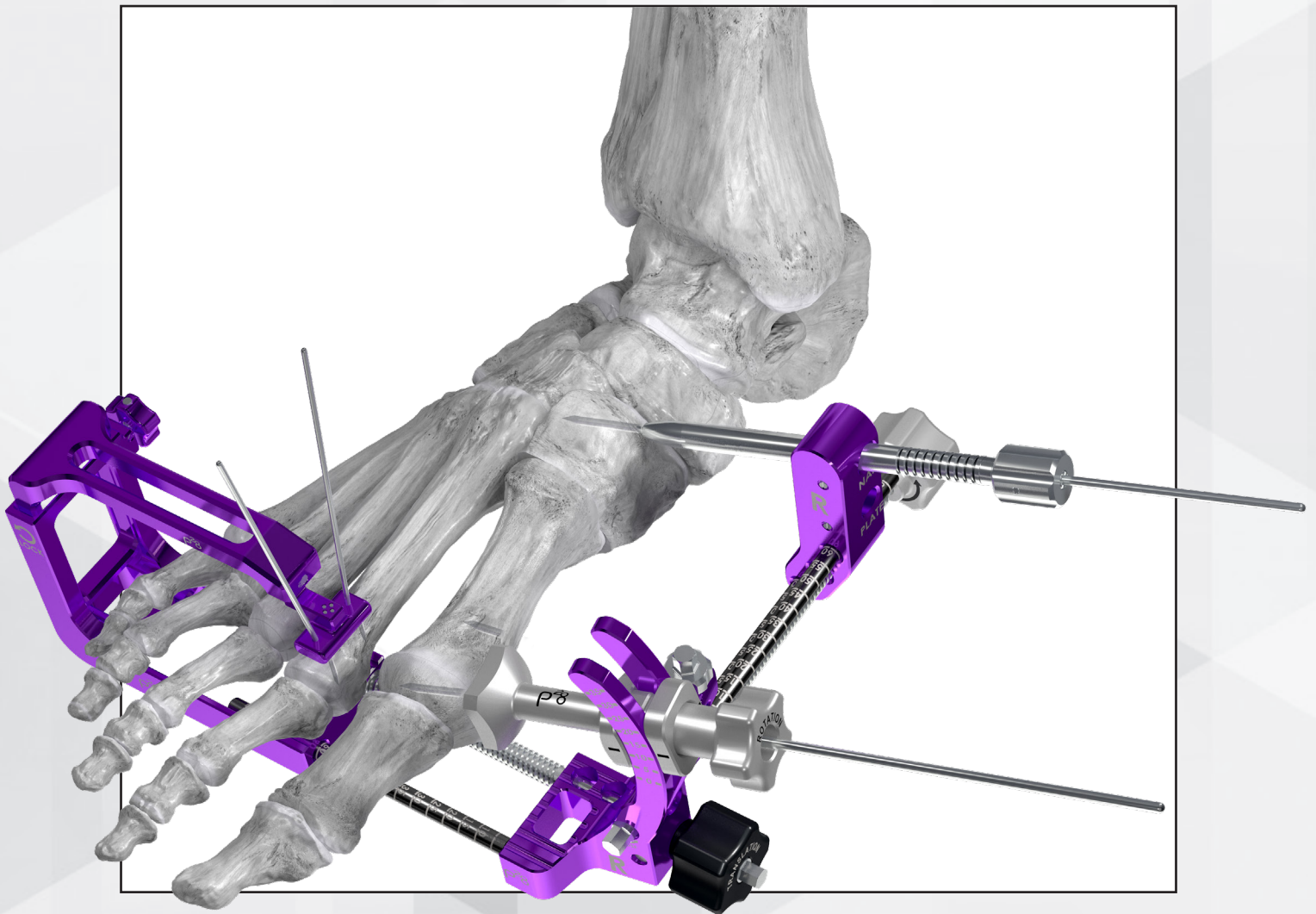


FEATURES OF THE BUN-YO-MATIC™ LAPIDUS CLAMP

- ▶ The Bun-Yo-Matic™ Lapidus Clamp is intended to remain on throughout fixation
- ▶ The clamp provides the ability to compress and distract based on surgeon need
- ▶ The Bun-Yo-Matic has multiple points of attachment to prevent slippage and grant the surgeon control of the foot
- ▶ The clamp is designed to avoid key visual landmarks to allow for interoperative fluoroscopy
- ▶ The Bun-Yo-Matic has the ability to de-rotate and translate the first metatarsal without the need for additional instruments or hands
- ▶ A foot wedge device is included to simulate weight bearing and allows the surgeon to load the foot hands-free

KEY FEATURES


Cut Guides

- ▶ Reproducible cuts can be achieved through the metatarsal and cuneiform cut guides
- ▶ Cuts are made based on the correction introduced by the clamp and designed to minimize resection and maintain length

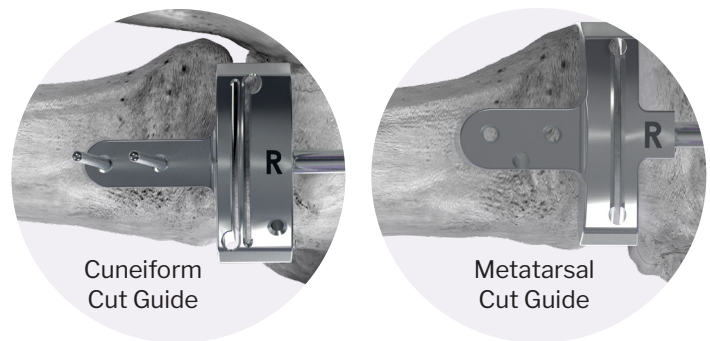


Intermetatarsal Angle Correction

- ▶ Translation of the first metatarsal can be performed to reduce the intermetatarsal angle
- ▶ The clamp provides up to 28° of intermetatarsal angle correction

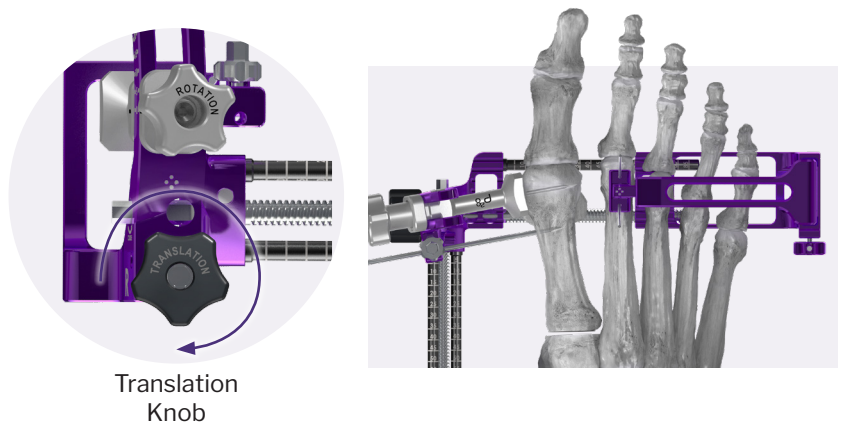

Metatarsal De-Rotation


- ▶ Grants up to 25° of first metatarsal de-rotation angle correction
- ▶ Rotates about the central axis of the metatarsal

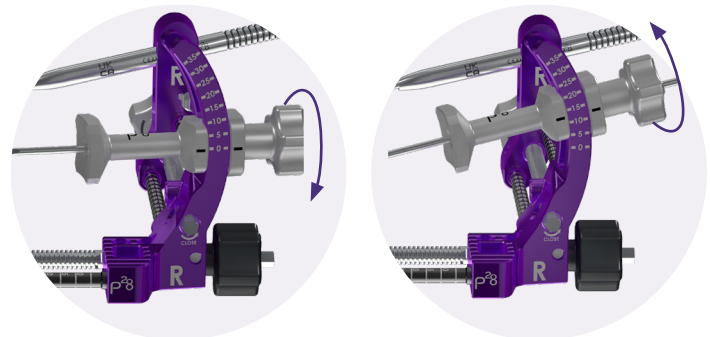


KEY FEATURES OF THE BUN-YO-MATIC™ LAPIDUS CLAMP

Joint Distraction and Closure

- Distraction for joint preparation and closure for alignment and reduction can be performed through a single device

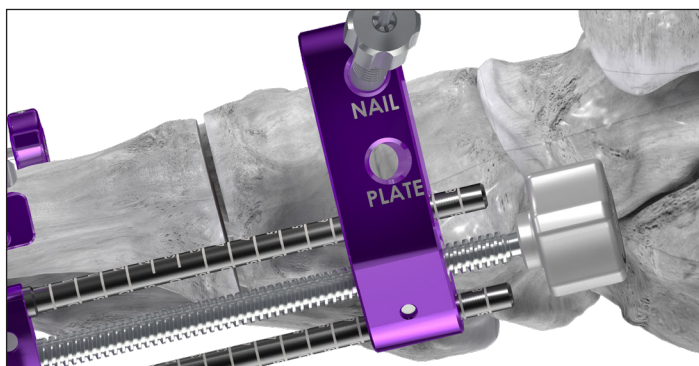



Before

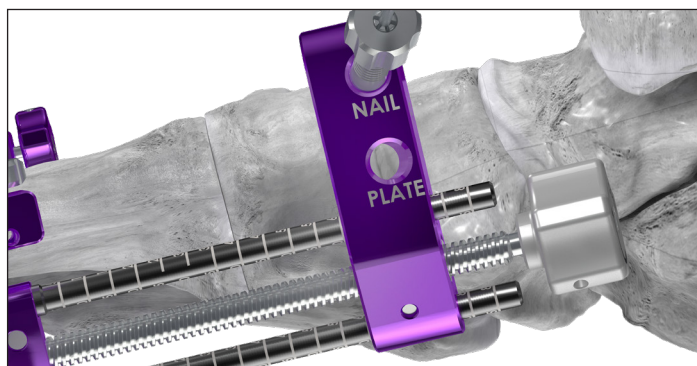

After

Multiple Fixation Options

- Select from Paragon 28's wide range of Lapidus Fixation options, including the Gorilla® Lapidus Arthrodesis System and the Phantom® Intramedullary Lapidus Nail

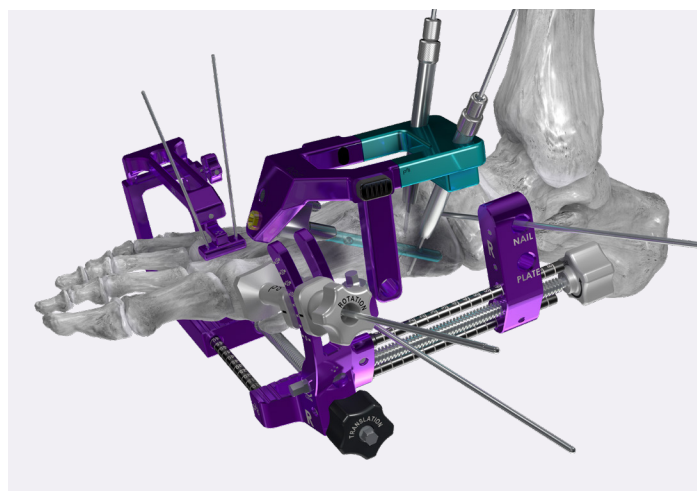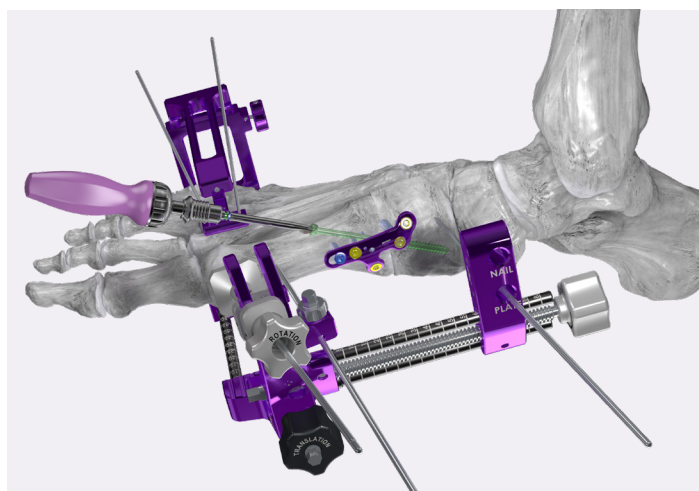

FEATURED INSTRUMENTATION:



K-wire, Ø1.40 mm x 10 cm



K-Wire, Olive, Threaded,
Ø1.80 mm x 4.7 cm



K-wire, Ø2.30 mm x 19 cm



Tissue Protector



Anti-Penetration Tool




K-wire Tube



Metatarsal Grip



Weight Bearing Wedge

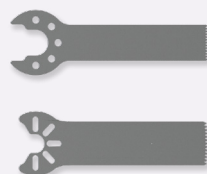


Saw Blades
(10 mm width)



Metatarsal Cut Guide



Cuneiform Cut Guide

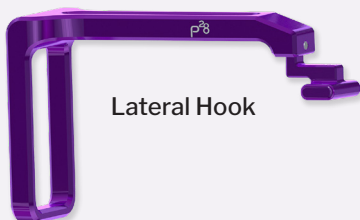


Lateral Hook



Metatarsal Paddle



Lateral Adjuster


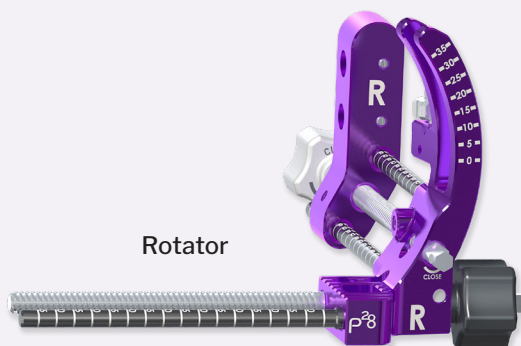






Rotator

TECHNIQUE SUMMARY

Clamp Attachment

Secure the clamp with K-wires

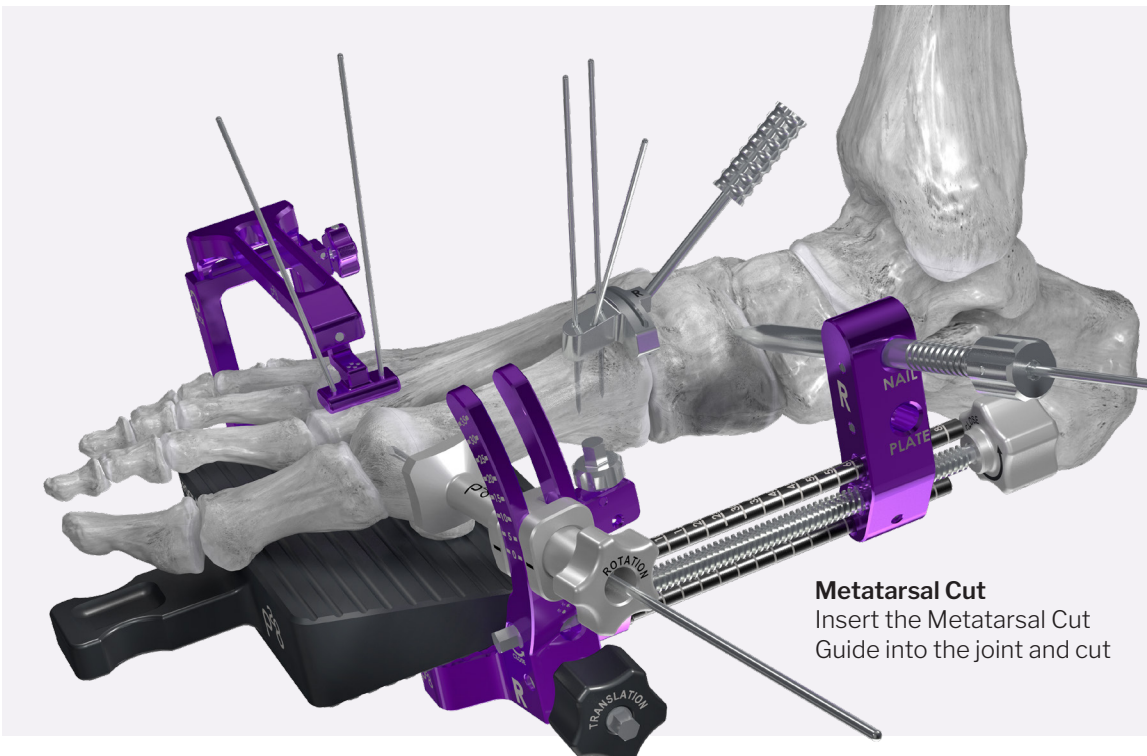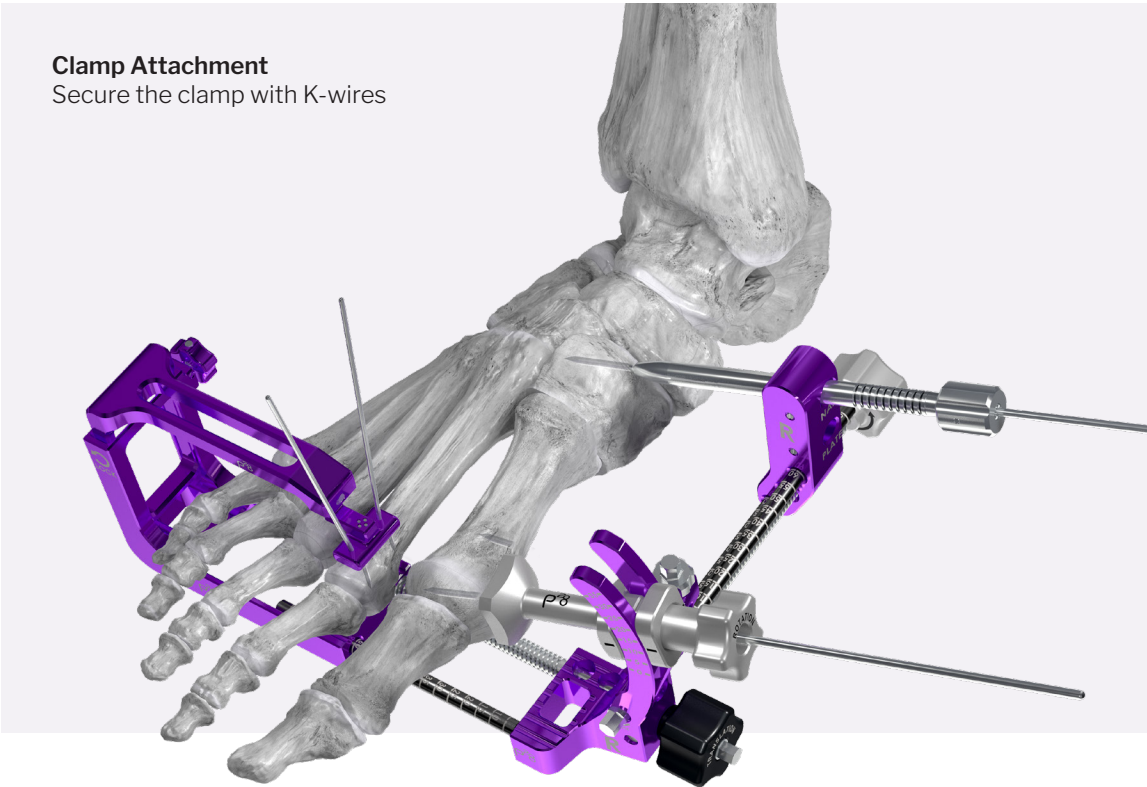

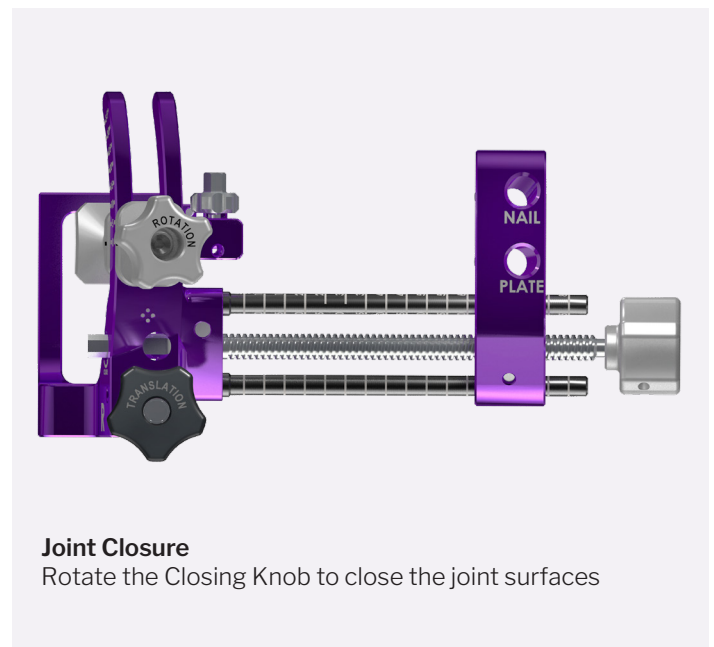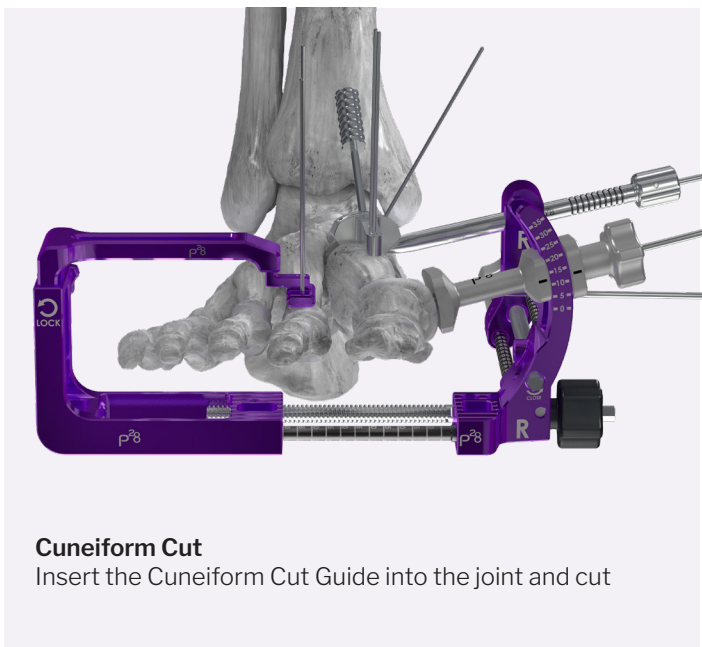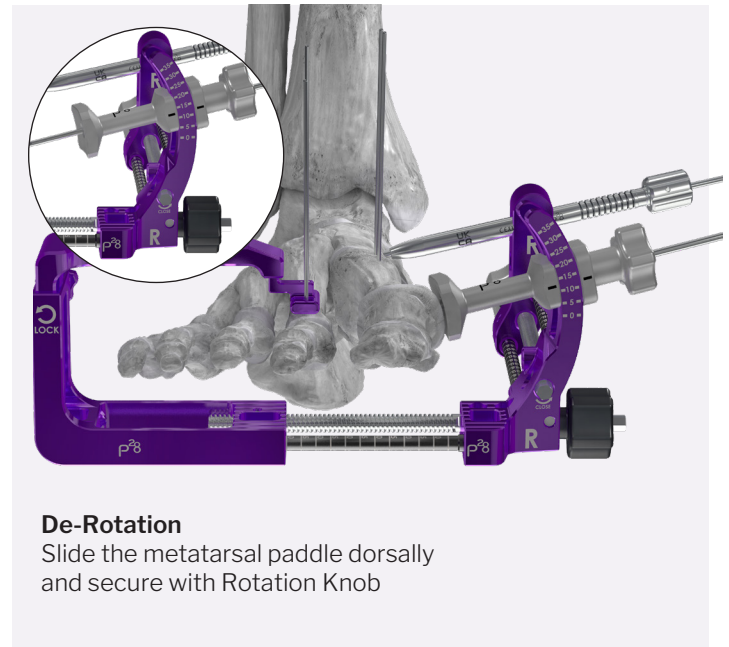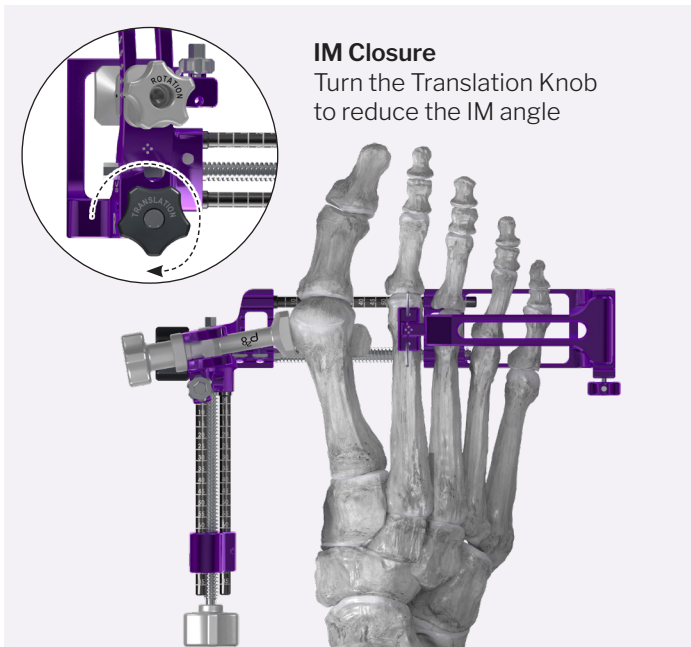




Bun-Yo-Matic™ Lapidus Clamp

BNYM-BRO-01 Rev A
2024-06-17

Paragon 28°, Inc. 
14445 Grasslands Dr.
Englewood, CO 80112 USA
(855) 786-2828

Paragon 28° Medical Devices Trading Limited
First Floor Block 7 Beckett Way
Park West Business Park
Dublin 12
D12 X884
Ireland
+353 (0) 1588 0350

™Trademarks and *Registered Marks of Paragon 28°, Inc.
© Copyright 2024 Paragon 28°, Inc. All rights reserved.
Patents: www.paragon28.com/patents

Exclusively foot & ankle ²⁸
Paragon®

www.Paragon28.com

For the contraindications, potential complications and adverse reactions, warnings and precautions associated with this device, please refer to the device specific instructions for use at <http://www.paragon28.com/ifus>