

Paragon²⁰[®]

Exclusively foot & ankle

BONE GRAFT[™] HARVEST SYSTEM

Bone Graft Harvest System PRODUCT INFORMATION



- Convenient instrumentation and straight forward technique to harvest cancellous autograft bone
- Available as a small, easy to use set compatible with the Gorilla[®] R3CON Plating System

Disposable Trepphine

Reusable Door

AO Quick Connect

Easy removal via twist knob to expose all graft for simple extraction from harvester

Secures morselized bone into housing until sufficient amount of graft is harvested


Depth markings in 10 mm increments

40 mm
30 mm
20 mm
10 mm

Morselizing tip allows consistency and volume of harvested bone to be controlled by surgeon

Bone Graft Harvest System SYSTEM CONTENTS

Bone Graft Harvest System

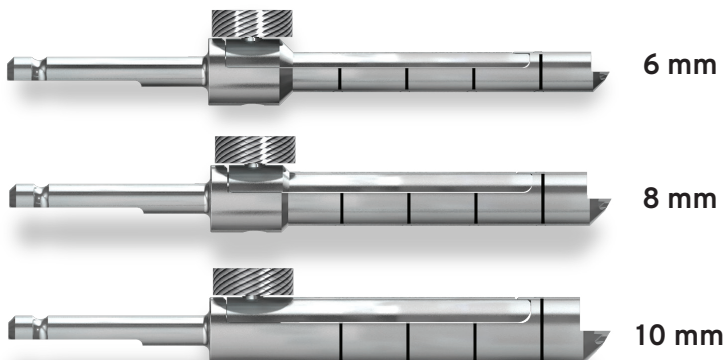

- Available in 6 mm, 8 mm and 10 mm diameters
- 2.3 mm K-Wire available for use if pilot hole is desired



K-Wire

Bone Graft Harvest System

ANCILLARY PRODUCTS


BONE GRAFT™

HARVEST SYSTEM

PRESERVE™ Bone Wedge System

Allograft Features:

- Aseptically processed
 - No gamma irradiation or bleach
- Density matched
- Available for Evans, Cotton, MTP, Lapidus, Subtalar, Calc-Cuboid and Hammertoe procedures

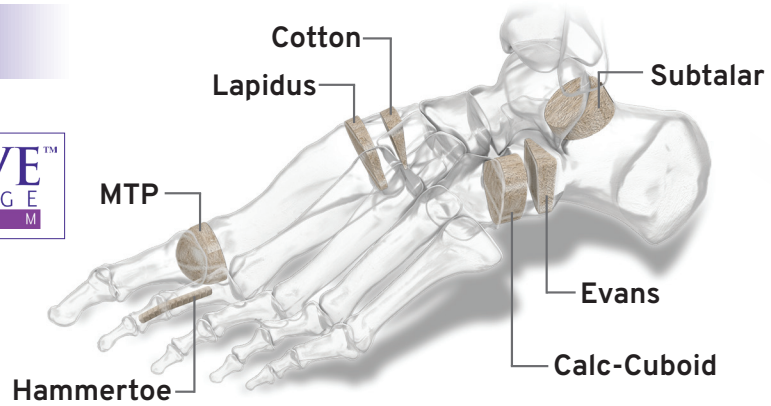

V92™ and V92 FC™ Cellular Bone Matrix

Allogenic™ bone matrix containing osteoprogenitor cells which aid in bone remodeling



- 92% Cell Viability Post Thaw
- Final handling properties yield a cohesive, wet sand consistency with V92™, and a moldable paste consistency with V92 FC™
- 150,000 Viable Cells/cc Available for Transplantation
 - Primitive, osteogenic cells only
- Novel Non-DMSO Cryoprotectant
 - No decanting required
- Optimized Bone Scaffold Component
 - Mineralized scaffold is a combination of cancellous and cortical bone
- Safe Tissue Processing
 - All graft portions from a single donor
 - Aseptically processed



OTHER PARAGON 28 BIOLOGIC SOLUTIONS



HARV-01 RevC

™ Trademarks and ® Registered Trademarks of Paragon 28, Inc.
 © Copyright 2020 Paragon 28, Inc. All rights reserved.
 Patents: www.paragon28.com/patents

Paragon 28, Inc. 14445 Grasslands Dr. Englewood, CO 80112 (855) 786-2828

Paragon 28 Medical Devices Trading Limited
 First Floor Block 7 Beckett Way
 Park West Business Park
 Dublin 12
 D12 X884
 Ireland
 +353 (0) 1588 0350



For additional information on Paragon 28® and its products please visit www.paragon28.com

