



## TITAN 3-D™ Wedge System

### PRODUCT INFORMATION

#### TITAN 3-D™ EVANS WEDGE

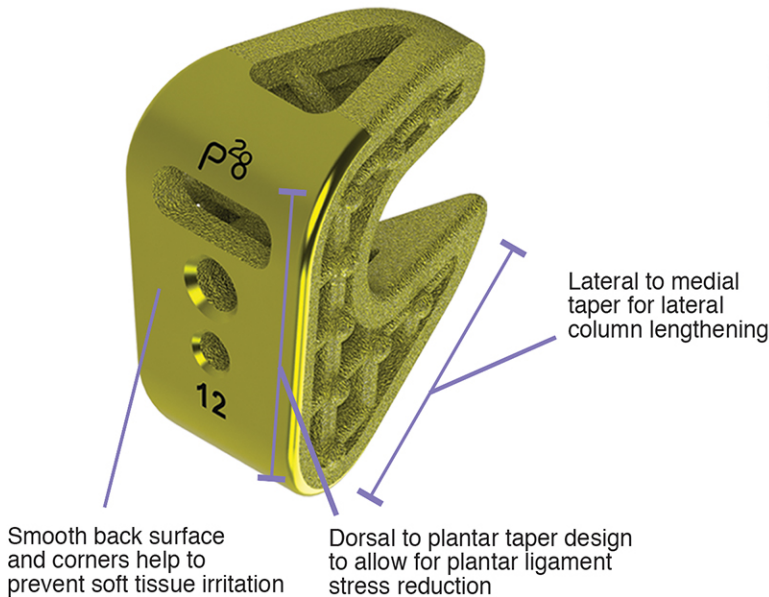


#### TITAN 3-D™ COTTON WEDGE

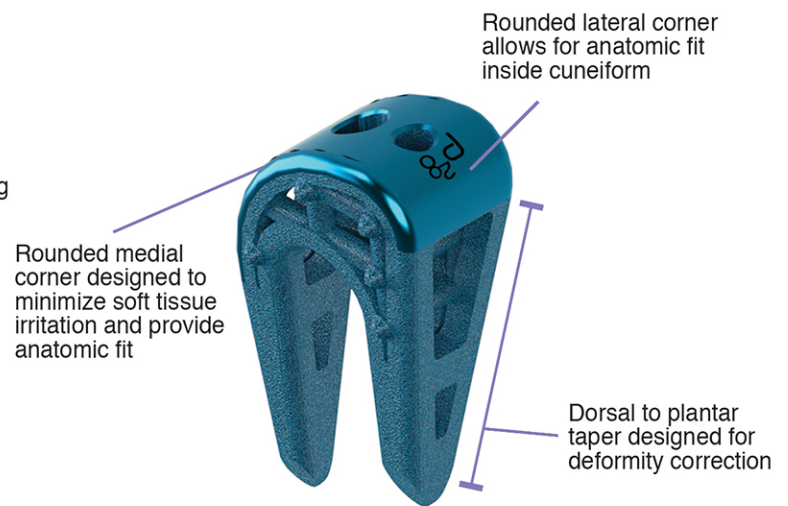


### FEATURES PRESENT IN BOTH IMPLANT FAMILIES

- Manufactured from medical grade titanium alloy (Ti-6Al-4V)
- Tapered nose helps to aid in implant insertion
- Open geometry with three-dimensional scaffold allows for blood entry, bone through growth and the incorporation of biologic products, if used
- Central opening allows for passage of screw across osteotomy to help increase stability of the construct
- Spikes on both sides of implant designed to help prevent expulsion from osteotomy site

### TITAN 3-D™ WEDGE SYSTEM - 12 Unique Titanium Wedges

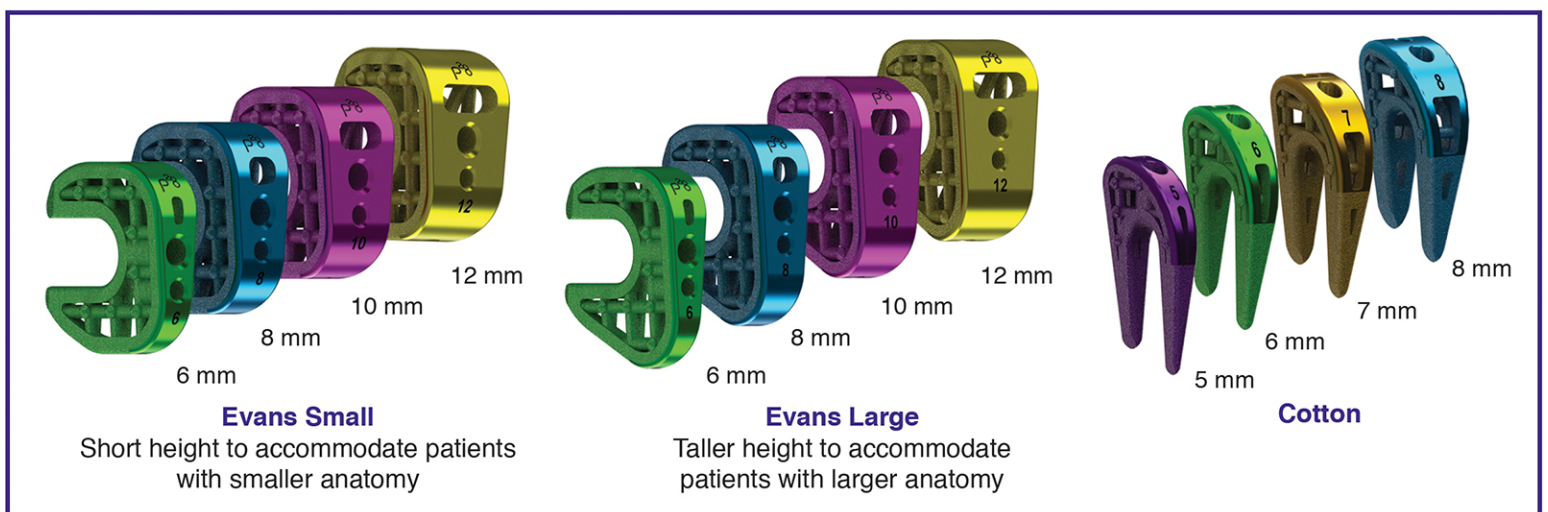



# TITAN 3-D™ Wedge System

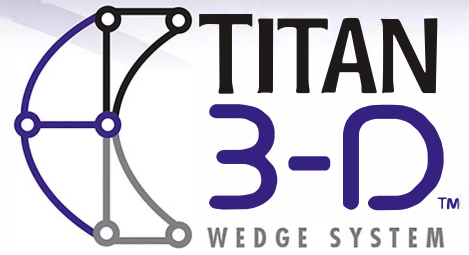

## PROCEDURE SPECIFIC INFORMATION

### PRODUCT SPECIFIC INSERTER (EVANS AND COTTON INSERTERS)

- Attaches to the back of the implant to facilitate accurate implant insertion and placement.

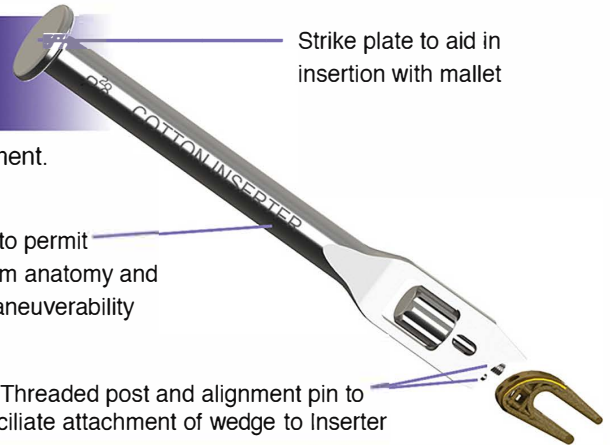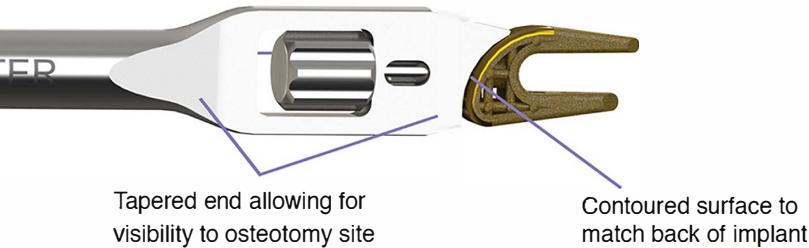


### PRECISION GUIDE SYSTEM (EVANS AND COTTON GUIDES)

- Designed to facilitate placement of screw through the wedge.



### RESECTION GUIDES (EVANS AND COTTON GUIDES)

- Available to aid in implant removal as necessary and designed to limit the removal of excessive patient bone

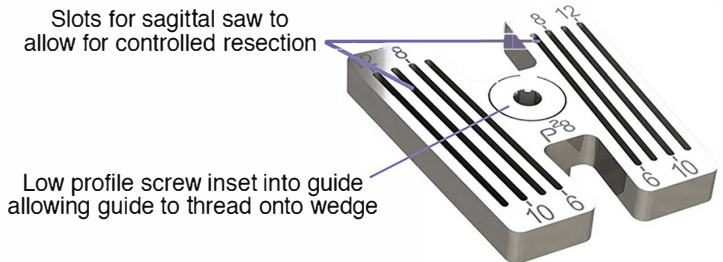

### ADDITIONAL FIXATION OPTIONS: (may be inserted through wedge)

#### 3.5 mm and 4.0 mm Mini-Monster® Fully Threaded Screws

3.5 mm Fully Threaded	4.0 mm Fully Threaded
10 mm - 50 mm (2.0 mm increments)	10 mm - 50 mm (2.0 mm increments) 55 mm - 60 mm (5.0 mm increments)

T3WS-01 RevC  
 ™ Trademarks and ® Registered Trademarks of Paragon 28, Inc.  
 © Copyright 2020 Paragon 28, Inc. All rights reserved.  
 Patents: [www.paragon28.com/index.php/patents](http://www.paragon28.com/index.php/patents)

For additional information on Paragon 28® and its products please visit: <http://www.paragon28.com>

Paragon 28, Inc.   
 14445 Grasslands Dr.   
 Englewood, CO 80112 USA   
 (855) 786-2828

Paragon 28 Medical Devices Trading Limited   
 First Floor Block 7 Beckett Way   
 Park West Business Park   
 Dublin 12   
 D12 X884   
 Ireland   
 +353 (0) 1588 0350

