

# JAWS<sup>®</sup>

NITINOL STAPLE SYSTEM

---

## BROCHURE

JAWS NITINOL STAPLE SYSTEM



Exclusively foot & ankle <sup>20</sup>  
Paragon<sup>®</sup>

## PRODUCT INFORMATION

The JAWS™ Nitinol Staple System uses superelastic nitinol and a simple insertion method to gain rigid compression across an osteotomy site. The staple is elastically deformed while stored in the inserter and immediately returns to its original shape after it is released from the inserter.


## IMPLANT FEATURES

- Sharp tooth geometry increases pull out resistance
- Lower surface profile than competitive designs to help minimize soft tissue irritation
- Staple achieves compression upon insertion; no heating required
- No plastic deformation required to achieve compression
- 15° angle staple option replicates the medial surface of the proximal phalanx
- Staple inserter helps eliminate the need to tamp
- Implant is biocompatible and doesn't need to be removed

## STAPLE SYSTEM OFFERING



8 mm  
Straight Staple

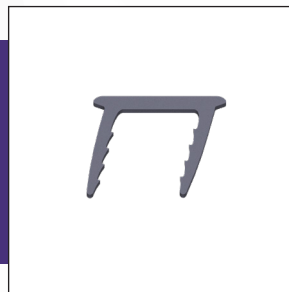

8 mm  
Angled Staple



10 mm  
Straight Staple



10 mm  
Angled Staple




12 mm  
Straight Staple




15 mm  
Straight Staple



18 mm  
Straight Staple

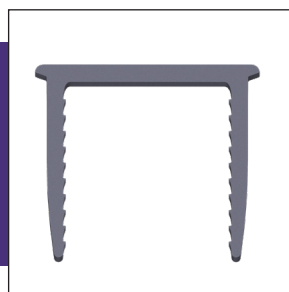

20 mm  
Straight Staple

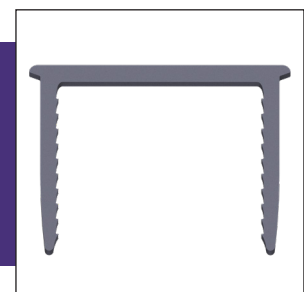

25 mm  
Straight Staple

## INSTRUMENTATION

The JAWS™ staple inserter is designed to allow the staple to be completely seated before it is released from the inserter. This allows for final placement of the staple before the staple fully compresses the osteotomy site.



Osteotome, 6 mm Straight



Osteotome, 12 mm Straight



Compressor

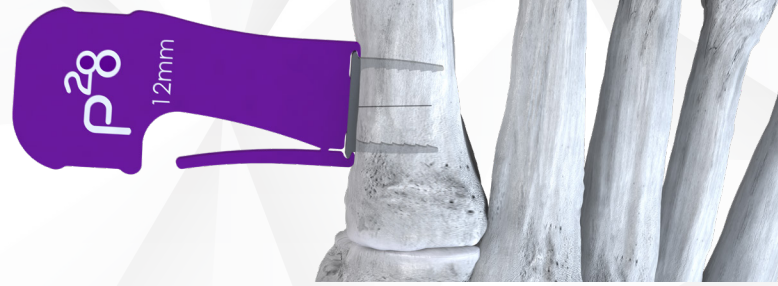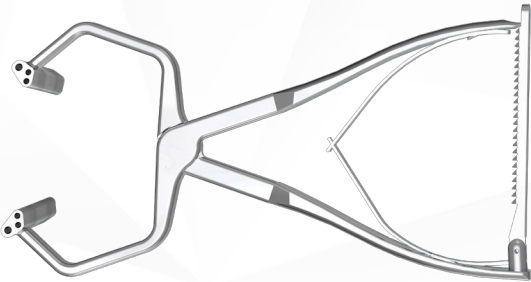

Osteotome, 19 mm Straight



Osteotome, 6 mm Curved



Osteotome, 12 mm Curved



Osteotome, 19 mm Curved




## STERILE PACKED INSTRUMENTATION

All instrumentation for the JAWS™ Nitinol Staple System comes in a simple, disposable, sterile kit. The staple comes pre-loaded on the inserter to help facilitate a quick and simple surgery.

8 mm, 10 mm or 12 mm Staple Kit

15 mm, 18 mm, 20 mm or 25 mm Staple Kit

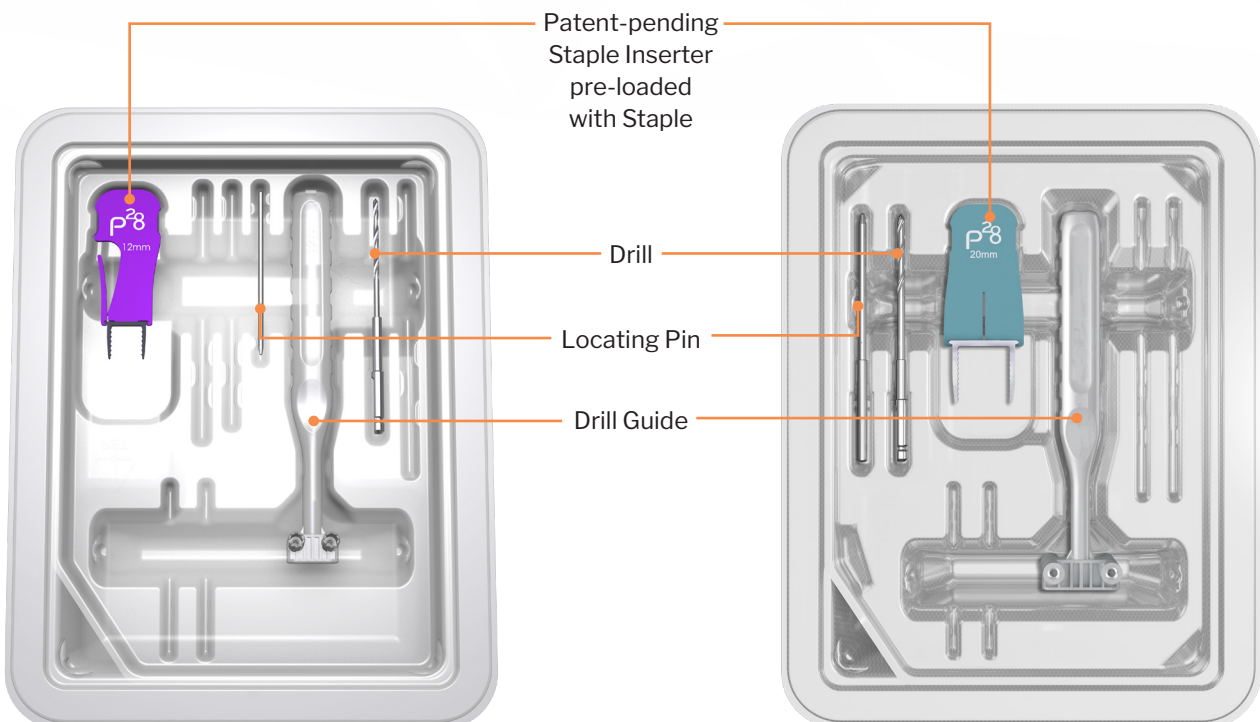

# JAWS<sup>®</sup>

## NITINOL STAPLE SYSTEM


---

PATENTED, DESIGNED & EXCLUSIVELY DISTRIBUTED BY

Exclusively foot & ankle <sup>28</sup>  
Paragon<sup>®</sup>

JNSS-01 Rev. E  
[2022/04/20]

<sup>™</sup>Trademarks and <sup>®</sup>Registered Marks of Paragon 28<sup>®</sup>, Inc.  
© Copyright 2022 Paragon 28<sup>®</sup>, Inc. All rights reserved.  
Patents: [www.paragon28.com/patents](http://www.paragon28.com/patents)

Paragon 28, Inc.   
14445 Grasslands Dr.  
Englewood, CO 80112 USA  
(855) 786-2828

Paragon 28 Medical Devices Trading Limited  
First Floor Block 7 Beckett Way  
Park West Business Park  
Dublin 12  
D12 X884

[www.Paragon28.com](http://www.Paragon28.com)